



VAN DER HEIJDEN
Labortechnik GmbH

c o o l s o l u t i o n s – y o u r a d v a n t a g e

Translation of the original operating instructions

Kühlmobil

311-B400-MK-FF-1-7-16-38-54

Artikel-Nr.: 1-100625

Edition: 01/11.2016

Foreword

CONGRATULATIONS

You have made a good choice. Thank you for your trust.

You have chosen a KÜHLMOBIL from Van Der Heijden Labortechnik GmbH.

These operating instructions is to familiarize yourself with the operation and capabilities of our KÜHLMOBIL.

Remember to follow all instructions and safety information.

Please abide by transportation, operation, maintenance, repair, storage and disposal to this manual.

By intended use of the KÜHLMOBIL, we offer full warranty and liability for our product!



Contents

1	Miscellaneous.....	6
1.1	Introduction.....	6
1.2	Instructions relating to Copyright and Trademark Rights.....	6
1.3	Instructions for the Operator	6
1.4	Refrigeration System Log Book	7
1.5	Service and Warranty	8
2	Safety	9
2.1	Miscellaneous	9
2.2	Correct Use	9
2.2.1	General references coolant secondary.....	10
2.2.2	General references coolant primary	12
2.3	Safety instructions in case of external water supply	14
2.4	Instructions on Signs and Symbols	14
2.5	Instruction Signs on the KÜHLMOBIL	15
2.6	Safety Instructions for Operating Personnel.....	16
2.6.1	Personal Protective Equipment	16
2.7	Safety Instructions for Operating the KÜHLMOBIL	16
2.7.1	Optical Check.....	16
2.7.2	Coolant Fluids	17
2.7.3	Leaking Testing.....	17
2.7.4	Mains Connection	17
2.7.5	Power supply with voltage switching	17
2.8	Safety Instructions for Maintenance	17
2.9	Notes on special types of hazard	18
2.9.1	Electric Power	18
2.9.2	Pumps	19
2.9.3	Coolant and other Chemical Substances	20
2.9.4	Hazards.....	20
2.9.5	Noise	22
3	Product Description	23
3.1	Explanation of Terminology	23
3.2	Overview	23
3.3	Functional Instructions.....	24
3.4	Accessories and Options	26
3.4.1	Floating Contact.....	26
3.4.2	Remote Control.....	26
3.4.3	Operating Hours Counter.....	26
3.4.4	Flow indicator	26
3.4.5	Flow control unit.....	26
3.4.6	External switch on.....	27
3.4.7	Flow monitor	27
3.4.8	Time Relay for Run-On	27
3.4.9	Overflow.....	28
3.4.10	Automatically re-fill.....	28
3.4.11	Bypass	28
3.4.12	Manometer	28
3.4.13	Dirt Filter	29

3.4.14	Fine Filter	29
3.4.15	Clean dirt filter	29
3.4.16	Clean fine filter	29
3.4.17	Interface RS232	29
3.4.18	Geodetic height difference	30
3.4.19	Tank drain	30
3.4.20	Automatic Switchover to Water Pipe in the Event of a Breakdown	30
3.5	Protective Equipment	31
4	Transportation, Storage and Commissioning	32
4.1	Transportation	32
4.2	Unpacking/Scope of Delivery	32
4.3	Put Up	33
4.4	Storage	34
4.5	Starting Up the KÜHLMOBIL	34
4.5.1	Bleeding the Circulating Pump	34
4.5.2	Cooling-Water Hoses	35
4.5.3	Primary	36
4.5.4	Floating Contact	36
4.5.5	External Switch-On Control	37
4.5.6	Remote Control	38
5	Touch-Pad Control Panel	39
5.1	Operating Features	40
6	Operating the Appliance	41
6.1	General	41
6.1.1	Operating Conditions	41
6.2	Switching 'ON'	41
6.3	Parametering and Switching 'ON'	42
6.4	Switching 'ON'	42
6.5	Continuous Operating Mode	43
6.6	Switching 'OFF'	44
6.7	Information Display	45
6.7.1	Parameter	45
6.7.2	Fault Memory	45
6.7.3	Appliance Data/Details	45
6.7.4	Display Settings	45
6.8	Fault Messages and Disruption Rectification	46
6.8.1	Schedule of possible Alarm States:	46
6.9	Disruption without Direct Message Reporting	48
6.10	Display Indications with Appropriate Options	49
6.10.1	Two pumps	49
6.10.2	External switching-on	49
6.10.3	The run-on time after external switching-on	49
6.10.4	Run-on time	49
6.10.5	Flow-Meter	49

7	Maintenance	50
7.1	Miscellaneous	50
7.2	Maintenance and Repair	50
7.3	Operating Status	51
7.4	Care and Visual Inspection	51
7.4.1	Visual Inspection	52
7.5	Servicing	52
7.5.1	General Instructions	52
7.5.2	Setting the Bypass Regulating Valve	53
7.5.3	Condenser (for Air Cooling)	53
5.4	Plate heat exchanger primary circuit	53
7.5.5	Tightness test according to EN 378	54
7.5.6	Coolant secondary circuit	54
8	Disposal	55
8.1	Environmental Protection	55
8.2	Coolant	55
8.3	Final Decommissioning	55
9	Following	56
9.1	Technical specifications	56
9.2	Primary water consumption	63
9.3	Components and spareparts	64

1 Miscellaneous

1.1 Introduction

These operating instructions are an essential aid for the successful and safe operation of the KÜHLMOBIL.

The operating instructions contain important directions for the safe, correct and economical operation of the KÜHLMOBIL.

Adhering to these instructions will help to avoid dangers, reduce repair costs and down-time and to increase the reliability and working life of the KÜHLMOBIL.

The operating instructions must be available at the KÜHLMOBIL at all times, and must be read and applied by every person engaged to work with/on the KÜHLMOBIL, e.g.:

- Operation, troubleshooting during work processing,
- Maintenance (servicing, care, repairs) and/or
- Transport.

1.2 Instructions relating to Copyright and Trademark Rights

These operating instructions are to be treated in confidence.

They should only be made accessible to authorised persons.

They should only be made available to third parties with the expressed, written permission of **Van der Heijden Labortechnik GmbH**.

All documents are protected in terms of the Copyright Act. Unless expressly allowed in writing, it is prohibited hand over or copy documents, even extracts thereof, or to sell or reveal their contents.

Any contravention is punishable and will incur damages.

We reserve all rights to exercise protection of industrial property rights.

1.3 Instructions for the Operator

The operating instructions form an essential part of the KÜHLMOBIL.

The operator must see to it that operating personnel are aware of these guidelines.

As the operating instructions manual will probably be subject to heavy use at the operating location, the operator is obliged:

- To keep the original safely,
- To ensure that a copy of the operating instructions is always available at the KÜHLMOBIL, and
- That every user reads the operating instructions carefully and follows and adheres to all instructions.

The operating instructions are to be supplemented by the operator to include operating instructions on the basis of existing, national regulations and accident prevention and environmental protection, including information on supervisory and reporting obligations relating to taking account of operational peculiarities, e.g. relating to work organisation and the personnel employed. In addition to the operating instructions and the regulations on accident prevention that are binding in the user's country and at the operating location, the recognised specialist regulations governing safety and correct work practices are to be observed.

The operator/user may not carry out any extensions/conversion to the KÜHLMOBIL, which could compromise safety, without the authorisation of **Van der Heijden Labortechnik GmbH**.

This applies in particular to the installation and adjustment of safety equipment.

Spare parts that come into use must meet the technical requirements stipulated by **Van der Heijden Labortechnik GmbH**.

This is always guaranteed for original spare parts.

Only use trained or instructed personnel. Clearly specify the responsibilities of personnel for operation, servicing and repairs!

Servicing of the KÜHLMOBIL may only be performed by specially trained **Van der Heijden Labortechnik GmbH**, personnel.

Van der Heijden Labortechnik GmbH reserves the right, in terms of product improvement, carry out technical modifications necessary to contribute to a proper function during the repair process.

1.4 Refrigeration System Log Book

From a refrigerant capacity of more than:

- 3.0kg → 134a or
- 1.3kg → R404a

the owner or operator must lead a logbook for this KÜHLMOBIL up to date according to EN378.

The system log (Logbook) must either be kept by the operator / owner and made available to the system or the information must be stored in a computer , with an print on the system !

In this case, the information must be available to the expert for maintenance or testing.

The Logbook includes :

- Details of all maintenance and repair work
- The quantity and nature of each time filled refrigerant
- Quantity and type of every time discharged refrigerant
- Changes and replacement of parts
- Results of all regular routine tests
- Longer downtimes
- Results of the annual audit of tightness

1.5 Service and Warranty

We endeavour to handle your questions and orders as quickly as possible.

We would ask you to indicate the identification data of your KÜHLMOBIL before every query. This data is entered on the cover sheet of these operating instructions.

Within the warranty period the **Van der Heijden Labortechnik GmbH** warranty for the KÜHLMOBIL covers damage proven to be as a result of design, material or manufacturing flaws.

This does not apply to any parts that are subject to wear and tear as a result of their function.

The warranty begins on the date of delivery by the manufacturer.

A complete warranty will only be honoured if:

- the KÜHLMOBIL is used in accordance with its intended purpose,
- the maintenance and service regulations are strictly adhered to.

The warranty is limited to the case of a complaint to our choice of the subsequent improvement or a free repair or a new delivery.

Defective parts will be repaired or replaced, provided that evidence present in the case of a disturbance or a defect of material or workmanship free of charge.

Any other compensation claims are excluded.

The spare-parts prices of **Van der Heijden Labortechnik GmbH** are applicable.

Subsequently we refer you to the warranty provisions as well as the general terms and conditions of **Van der Heijden Labortechnik GmbH**.

Please direct any queries to **Van der Heijden Labortechnik GmbH** You will find the address on the title page.

2 Safety

2.1 Miscellaneous

The KÜHLMOBIL is constructed in accordance with state-of-the-art technology and recognised safety regulations.

Nevertheless use of the unit can involve dangers for the user or third parties as well as adverse effects on the KÜHLMOBIL itself or other material assets if it:

- Is operated by untrained or non-familiarised personnel,
- Is not used correctly,
- Is maintained or serviced incorrectly.

2.2 Correct Use

The KÜHLMOBIL by **Van der Heijden Labortechnik GmbH** is used exclusively for cooling externally connected devices and machines, which need to be cooled at a constant supply temperature.

The cooling medium between the KÜHLMOBIL and the devices being cooled is “HKF 10.1.” or distilled water with sodium bicarbonate (2.2.1).

Deionised water or distilled water may only be used in pure form as a cooling medium when all components in the cooling circuit made of stainless steel or plastic.

The coolant HKF 15 is to be used, with KÜHLMOBILen which were conceived for outdoor installation and the ambient temperature can amount to -15° C at a stand still.

Inlet temperatures <10° C needed also an anti-freeze protection with HKF 15.

Are at a standstill of the KÜHLMOBIL ambient temperatures below -15 ° C expected, please ask at **Van der Heijden Labortechnik GmbH** for the accordingly cooling fluid.

The KÜHLMOBIL may only be connected to devices that are water-cooled and correspond to the parameters of the KÜHLMOBIL.

The thermal load of the unit to be cooled must not exceed the cooling output of the cooler.

The through-flow resistance of the unit to be cooled must be approx. 15% below the maximum indicated pump pressure.

Only those devices to be cooled, which are authorised for the maximum water-pump pressure, may be connected.

A different use or use that goes beyond this description, e.g. as a refrigeration system for food, is considered incorrect.

The operator of the KÜHLMOBIL shall bear sole liability for any damage resulting from incorrect use.

This applies similarly to unauthorised modifications to the KÜHLMOBIL. Components/spare parts that are not in fault-free condition are to be replaced immediately.

Only use original spare parts. Under certain circumstances correct function cannot be guaranteed if non-original spare parts are used!

Correct use also include following the following instructions:

- on safety,
- on operation and
- on maintenance and
- for servicing

which are described in these operating instructions.

Attention! If there is a defect in the cooling fluid circuit there is the risk that oil from the refrigeration cycle (compressor) of chiller can get into the coolant!

Observe all laws and regulations of the water-supply company in force at the operation site!

2.2.1 General references coolant secondary

All KÜHLMOBILs, which are equipped with transfer pumps may only be operated with the cooling liquid "HKF 10.1".

This applies for temperatures in the secondary water flow $\geq +10^{\circ}\text{C}$.

The coolant can be obtained from the Fa. **Van der Heijden Labortechnik GmbH**. At lower secondary flow temperatures asking the manufacturer.

As an alternative to "HKF10.1" can also be a mixture of distilled water and sodium hydrogen carbonate are used. 1.68 grams of NaHCO_3 is to mix per litre of distilled water.

Attention! The mixing ratio must be observed, as a PH in this ratio of the amounts - should set the value of approximately 8.0.

This target PH - value must be monitored and checked every six months.

The neutral pH - value and lime-free and residue-free water can settle to the mechanical seal a feed pump nothing prevents leakage over a prolonged period.

Attention! For magnetically coupled pumps This is not true.

All KÜHLMOBILs, which are equipped with demineraliser may only be operated with demineralized water. This applies for temperatures in the secondary water flow $\geq +10^{\circ}\text{C}$ to.

All water-carrying components in the secondary cooling circuit of KÜHLMOBIL are made of plastic or stainless steel. The same applies to the operation of the automobile with distilled water cooling.

The neutral pH - value and lime-free and residue-free water can settle to the mechanical seal a feed pump nothing prevents leakage over a prolonged period.

When using normal drinking water or well water with blended Demineralized water voids any warranty on the pump for leaks if build-up on the mechanical seal are determined by using the above liquid.

In sealless (magnetically coupled) pump this point does not apply.

A defective mechanical seal must therefore always be sent in for Van der Heijden Labortechnik GmbH, so that appropriate examinations can be carried out.

Likewise the cleanliness of each medium must be guaranteed. The algae retardants "Thermo Clean" may be used as an additive for water.

The dosage is given in the data sheet and must be strictly observed.

(Thermo Clean can be ordered from the Van der Heijden Labortechnik GmbH.)

Light through all water-carrying parts must be avoided, as this promotes the growth of algae in the medium. By algae in the water, the pumps have a mechanical seal leaking.

The mechanical seal of the pump is a consumable item. Therefore, deposits by example hard water, improper additives, growth of algae and bacteria, etc. to be avoided.

Note: deposits inevitably lead to leaks and malfunction of the pump.

Leaks at the pump can also arise when commissioning has not been carried out according to regulations.

Commissioning is mandatory to carry out according to the instructions as it may result in a short-term dry running of the pump otherwise.

By running dry are scoring form at the mechanical seal. The result is a leaky pump.

The use of the following cooling media (also mixed with water) is deemed not determination in accordance with:

- Solvents of any kind,
- Flammable liquids,
- Explosive liquids,
- Caustic liquids

2.2.2 General references coolant primary

(for all water cooled models)

With water-cooled KÜHLMOBILs from any kind, which will be connected on a drinking water network, house cooling system or something else, the data sheets for the composition of the necessary cooling water have to be taken from the manual

Despite of high-quality V4A heat exchangers an attack on material or on larger connections can be in disregard. A primary water filter has to be installed before the water inlet.

The filter size should be 90 µm at least. Rust particles of steel lines can also destroy V4A-plates by corrosive behavior. Guarantees are not possible in these cases.

The cause of pitting or something like that results high costs. This applies first of all to KÜHLMOBILs with active cooling.

Fouling

It has to be ensured that the DIN guidelines for drinking water and heating water, Vd-TÜV guidelines for ADFW as well as WTT guidelines have to be observed (see diagram).

Many different factors can influence fouling. These are, for example: velocity, temperature, turbulence, distribution, water quality.

The mediums have to be moved at the highest possible mass flows. In the event of excessively low mass flows (part load) the turbulence in the heat exchanger can decrease and the fouling tendency increases. Lime deposits in the heat exchanger can be possible at temperatures above 60° C according to the water quality.

Turbulent flow and lower temperatures reduce the risk of calcification. During switching-off the KÜHLMOBIL it has to be ensured that the secondary side will be switched-off first and then the primary side.

During switching-on the primary side has to be switched-on first and after that the secondary side. In that way overheating of the heat exchanger will be avoided.

Cleaning

If it should be expected a plaque formation based on the water quality (for example high degree of hardness or strong pollution), a cleaning has to be performed in regular distances.


There is the possibility of cleaning by rinse. In that case the heat exchanger has to be rinsed with suitable cleaning solvent against the normal way of flow.

If chemicals will be used for cleaning, please note, that there is no indigestibility against stainless steel, copper or nickel. Non-observance can cause the destruction of the heat exchanger.

Corrosion resistance of soldered plate heat exchangers compared with water ingredients.

The brazed plate heat exchanger consists of embossed stainless steel plates 1.4404 or SA240 316L. It is thus to be considered the reaction of corrosion of stainless steel and solder copper or nickel.

Water constituent + parameters	Unit	Heat exchanger, copper brazed	Heat exchanger, nickel brazed
pH-value		7 – 9 (considering SI Index)	6 - 10
Saturation-Index SI (delta pH-value)		-0,2 < 0 < +0,2	No specification
Total hardness	°dH	5-7	5-7
Conductivity	µS/cm	50...250	No specification
Abfilterbare Stoffe	mg/l	<20	<20
Chlorides	mg/l	<25	
Free Chlorine	mg/l	<0,3	<0,3
Hydrogen sulphide (H ₂ S)	mg/l	<0,05	No specification
Ammonia (NH ₃ /NH ₄)	mg/l	<2	No specification
Sulphates	mg/l	<50	<150
Hydrogen carbonate	mg/l	<200	<200
Hydrogen carbonate / Sulphates	mg/l	>1,0	No specification
Sulphide	mg/l	<1	<5
Nitrate	mg/l	<50	No specification
Nitrite	mg/l	<0,1	No specification
Iron	mg/l	<0,2	No specification
Manganese	mg/l	<0,1	No specification
Free aggressive carbonic acid	mg/l	<15	No specification

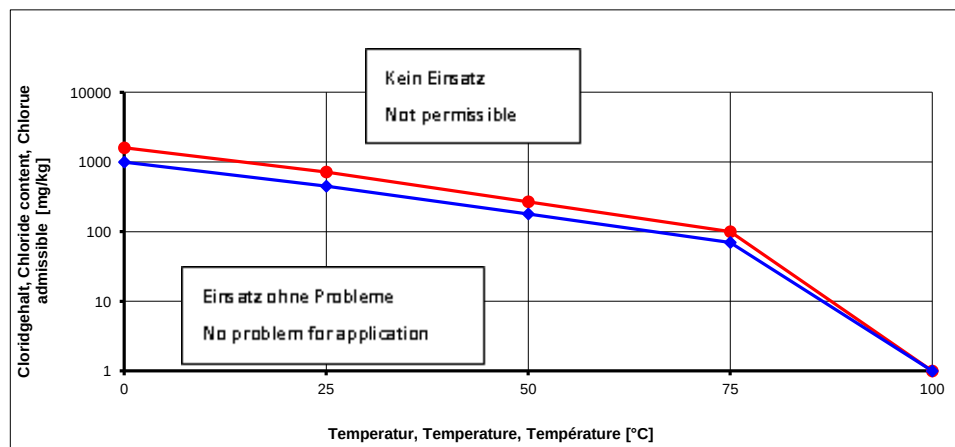

2.3 Safety instructions in case of external water supply

Due to the design of the as a water-cooled or water - water – version the KÜHLMOBIL is supplied with cooling fluid of the domestic water system or municipal water on the primary circuit. The same applies to the options "emergency cooling" and "automatic refill". Both options must be connected to their function to an external water circuit.



Since the external water supplies are not self locking, it can lead to leakage in the event of an uncontrolled water leakage.

For this reason, the corresponding area on the KÜHLMOBIL, at the hose connections and at the device which have to be cooled must be monitored with moisture sensors.

An additional floor drain is recommended.

By a leakage to avoid further damages, it must be ensured that all water supplies and the KÜHLMOBIL shut off immediately and be protected against restart!

2.4 Instructions on Signs and Symbols

In these operating instructions the following designations and signs are used for safety instructions and particularly important details:



Danger!

Draws attention to an immediately threatening danger, which will lead to the most serious personal injuries or to death if the indicated instruction is not followed exactly.



Warning!

Draws attention to a possibly dangerous situation, which could lead to the most serious personal injuries or to death if the indicated instruction is not followed exactly.



Caution!

Draws attention to a possibly dangerous situation or unsafe and dangerous methods of procedure, which could lead to personal injury or to material damage to the KÜHLMOBIL or its environment.

All of these WARNINGS must be precisely adhered too!

This refers to particular circumstances, adherence to which guarantees safe, correct and efficient handling of the KÜHLMOBIL. All instructions should be followed in the interests of correct use of the KÜHLMOBIL. Pass on all instructions to other users as well!



Note!

This refers to particular circumstances, adherence to which guarantees safe, correct and efficient handling of the KÜHLMOBIL. All instructions should be followed in the interests of correct use of the KÜHLMOBIL. Pass on all instructions to other users as well!

- The bullet identifies work and/or operating steps. The steps are to be performed in sequence from top to bottom!
- The mirror hyphen identifies enumerations.

Instructions and symbols such as warning signs, warning stickers, brief instructions, component identifications etc., that are attached directly to the KÜHLMOBIL, absolutely must be adhered to.

Instructions and symbols attached directly to the KÜHLMOBIL must not be removed, and must be maintained in fully legible condition!

2.5 Instruction Signs on the KÜHLMOBIL

Following instruction signs are on the chiller:

BREIF INSTRUCTION

ATTENTION! Follow before starting the device:

1. The first connection with the adjustment of the operating voltage must be performed by a qualified electrician.
2. It is essential to bleed the water pump
3. Fill the KÜHLMOBIL container with coolant up 2,5 cm to the rim.
4. Attach a piece of hose to the water supply and hold it in a bucket.
5. Open the advance and wait till the coolant flows into the bucket bubble-free.
6. Now the pump bleeds the system automatically, adding coolant may be necessary.
7. Please follow the operating instructions!

This sticker goup to **size 08** inside lift-up cover next to the KÜHLMOBIL rating plate.

BREIF INSTRUCTION

ATTENTION! Follow before starting the device:

1. The first connection with the adjustment of the operating voltage must be performed by a qualified electrician.
2. It is essential to bleed the water pump
3. Fill the KÜHLMOBIL container with coolant up 7.0 cm to the rim.
4. Attach a piece of hose to the water supply and hold it in a bucket.
5. Open the advance and wait till the coolant flows into the bucket bubble-free.
6. Now the pump bleeds the system automatically, adding coolant may be necessary.
7. Please follow the operating instructions!

This sticker is **from size 1** attached backside the chiller next to the KÜHLMOBIL rating plate.



This sticker is attached above the mains cable outlet at the rear of the KÜHLMOBIL.

**externe
Einschaltung**

These stickers with different text content naming the respective component are attached to all fittings, connections and cables at the rear of the KÜHLMOBIL.



The CE symbol is also attached at the rating plate of the KÜHLMOBIL.

Rating Plate

The rating plate is attached to the inside of the KÜHLMOBILs; fold-up cover.

You will find a picture of the rating plate in the technical specifications in the appendix of these operating instructions.

2.6 Safety Instructions for Operating Personnel

The KÜHLMOBIL may only be used in a technically fault-free condition as well as in accordance with its intended purpose and in a safety- and danger-conscious manner in accordance with these operating instructions! Any faults, especially those that can negatively influence safety, must be rectified immediately!

Anyone who is engaged with setting-up, starting-up, maintenance or troubleshooting associated with the KÜHLMOBIL must have completely read and understood these operating instructions – especially the *Safety* chapter *before commencing work*.

It is to late during operational use.

This is particularly true for personnel who only occasionally operate the KÜHLMOBIL.

The operating instructions must be kept ready to hand at the KÜHLMOBIL at all times.

No liability will be accepted for damages and accidents resulting from failure to follow the operating instructions. The relevant accident-prevention regulations as well as the other generally recognised safety-related and work-related medical rules must be observed.

Should malfunctions or safety-related changes become apparent in the behaviour of the KÜHLMOBIL, it should be shut down immediately and the event should be reported to the responsible office/person!

Work on the KÜHLMOBIL may only be performed by reliable and trained personnel.

Only use trained or familiarised personnel!

2.6.1 Personal Protective Equipment

Personal protective equipment should be for each person who is involved in maintenance, repair and recovery, ready made and used:

- For all refrigerants and independent of the characteristics of the refrigerant:
- ***protective gloves, safety shoes (S3) and eye protection***

2.7 Safety Instructions for Operating the KÜHLMOBIL

For all work, which relates to

- operation,
- adjustment of the KÜHLMOBIL and its safety equipment and
- maintenance,

Switching on and off must be done in accordance with these operating instructions and the maintenance instructions must be adhered to!

The KÜHLMOBIL may only be operated when all protective and safety equipment, e.g. detachable protective equipment, housing locks, etc. are present and functioning!

Never operate units with damaged mains power cables.

2.7.1 Optical Check

At least every six months must KÜHLMOBIL be checked for visible damage. Changes (including the operating performance) immediately report to the competent shift supervisor or the superintendent.

2.7.2 Coolant Fluids

The coolant quality must be checked at regular intervals (min. Once a month) to avoid contamination.

Chemical or biological contaminants must be removed immediately to damage to the KÜHLMOBIL, especially at the pump to avoid.

2.7.3 Leaking Testing

The refrigeration cycle with appropriate capacity is once a year to check means of leak test according to EN 378-2. (see also under item 1.4 and 7.5.5)

The tightness test is responsible for an authorized service technician or qualified person in compliance with all safety regulations and must be recorded in the logbook.

2.7.4 Mains Connection

The power plug serves as safe disconnecting device from the mains and must be in-mer freely accessible.

The unit never operate units with damaged mains power cables.

2.7.5 Power supply with voltage switching

Before the first terminal of KÜHLMOBIL to the supply voltage of the chip voltage selector switch must be set by authorized technicians for installation.

Factory the KÜHLMOBIL for an operating voltage of 230VAC (switch position 0) is set.

Subsequent adjustment at a later time, for example, by the operator is under-says. Only in this way damage to the KÜHLMOBIL can be avoided.

Again 2.7.4 !!

2.8 Safety Instructions for Maintenance

Inform the operator of the device to be cooled before the start of the implementation of special and maintenance work.

The prescribed or indicated in this manual deadlines for recurring inspections / compliance inspections.

At inspection, maintenance and repair of KÜHLMOBIL strictly the instructions for maintenance work in Chapter 7, Maintenance note s and the safety devices!

To carry out maintenance work, the personal protective equipment and work tools and workshop equipment is absolutely necessary.

The KÜHLMOBIL must be switched off for maintenance, repair and repair work and secured against unexpected reconnection by:

- The KÜHLMOBIL is switched to the standby mode with the button  and then
- The power cable of KÜHLMOBIL s (eg by pulling plug from the socket) is disconnected from the mains.

- Cool down hot components after turning off the KÜHLMOBIL.

Work on electrical equipment may only be taken before-only by qualified electricians!

For maintenance and repair work to loosened screws connections then always re-tighten!

If the dismantling of safety devices during maintenance or servicing work required must be carried out immediately after the completion of maintenance and repair work, the assembly and testing of the safety devices!

The KÜHLMOBIL, and especially connections and fittings are to release at the beginning of the maintenance / repair / maintenance of all dirt and residue, such.

As dust, supplies or maintenance products.

Only use mild detergents, water-based. Observe the manufacturer. No organic solvent used as fire and explosion hazard!

2.9 Notes on special types of hazard

2.9.1 Electric Power

Only use original fuses with specified current on the mains side!

Never perform work on voltage-carrying components.

When carrying out repairs, ensure that design features are never altered in a way that might compromise safety. In particular, creepage and air ways as well as gaps must not be reduced due to insulation.

Fault-free earthing of the electrical system must be guaranteed by means of a protective conductor system.

In the event of interruptions in the electrical power supply, immediately disconnect the KÜHLMOBIL from the mains by:



The perfect grounding of the electrical system must be ensured by a protective grounding system.



An automatic disconnection of the current supply by a fault current protection device must be present, in order to prevent a condition danger-bringing by a contact voltage.

In the event of interruptions in the electrical power supply, immediately disconnect the KÜHLMOBIL from the mains by:

- Switching the KÜHLMOBIL to standby mode by means of the  key, and then
- Unplugging the mains cable of the KÜHLMOBIL from the socket.

2.9.2 Pumps

References to modifications, repairs or changes:

References at the pump consider and keep it readable, e.g. direction of rotation arrow, marking of the coolant connection.

Don't move parts during the enterprise depart for contact protection for hot and cold.

If necessary use protection equipment.

Work only on the pump when it is off.

At all assembling and maintenance work scoled the motor zero potential and protect it for against restarting.

Install the safety devices again correctly on the pump after all work.



Residual Risk:

Draw in from long, open hair at the protective covers of the pump is possible.

- Carry a hairnet!

Injuries caused by expel objects from the opening of the fan motor cover, which are introduced to the opening.

- Do not introduce any objects!

Burns or scaldings by hot surfaces or hot coolant.

- Do not touch!
- Weare safety gloves!

Injuries caused by cooling liquid contract out of defect sealing.

- Take the pump out of operation!
- Repair the pump!



Special risks:

Magnetic drive

The strong magnetic field within the range of the magnetic couplings or with single magnets can lead to following dangers:

- Mortal danger for persons with cardiac pacemaker
- Damages at magnetic data media (document of identification with magnetic stripes, credit and cheque cards), electrical, electronic, fine-mechanical devices, construction units or instruments (e.g. mechanical, digital clocks, pocket calculator, non removable disks)
- Uncontrolled tightening of magnetic parts (e.g. tools, screws)



NOTE!

- Keep a safety distance of minimum 150 mm to open magnets or magnetic couplings with magnet-sensitive objects.
- For plug-in units which are not installed into the drive unit and not completed drive units keep a safety distance of minimum 150mm for magnet-sensitiv parts and cardiac pacemaker



NOTE!

Distance to the completely installed pump:

The magnetic fields of the magnetic couplings by installed pumps become completely shielded by the enclosing construction units. By magnetic coupling goes neither in the stop nor in the enterprise a danger out.



DANGER!

Strong magnetic field within the range of the magnetic coupling, with single magnets, plug-in units or pump heads.

Mortal danger and damage to property caused by the magnetic field!

- Guarantee that personnel with an cardiac pacemaker does not accomplish work on the pump.
- Protect the workstation, if necessary block the workstation off:
 - Guarantee that personnel with an cardiac pacemaker have > 1 m safety distance.
 - Guarantee that no magnetizable metal parts can be tightened by the magnetic coupling.
 - Guarantee that parts of the magnetic coupling cannot be tightened by magnetizable metal parts
- With magnet-sensitive articles keep a safety distance > 150 mm to the magnetic coupling or pump head.

2.9.3 Coolant and other Chemical Substances

When handling chemical substances, observe and adhere to the relevant regulations and safety specifications sheets of the manufacturer of these substances relating to storage, handling, use and disposal!

Eating, drinking, smoking and keeping food in rooms in which chemicals are located is prohibited!

Never keep hazardous materials in food containers or receptacles.

Always use and clearly label containers authorised for the respective material.

With regard to the coolant, please adhere to the safety specifications sheet accompanying the device/operating instructions.

2.9.4 Hazards

Become vaporizer or air coolers installed in the immediate vicinity of heat sources, to take effective measures to prevent these vaporizer or air coolers are exposed to too much heat, which leads to high pressures.

Condenser and liquid receiver must not be positioned in close proximity to heat sources.

Threats by the ruling in refrigeration temperature conditions are mainly due to the simultaneous

occurrence of liquid and vapor phases.

In addition, the state of the refrigerant and the stress on the various parts not only on the operations and functions within the system depends on, but also by external factors.

The following hazards are noteworthy:

- *Hazards by direct influence of extreme temperatures such as:*
 - Material embrittlement at low temperatures
 - Freezing of trapped liquid
 - Thermal Stress
 - Changes in volume caused by temperature changes
 - Personal injury, caused by low temperatures
 - Accessibility of hot surfaces

- *Hazards caused by increased pressure e.g.:*
 - Increase in the liquefying pressure caused by insufficient cooling
 - Increase the vapor pressure by heating from the outside
 - Expansion of liquid refrigerant due to increase in the outside temperature
 - Fire

- *Hazards by direct influence of the liquid phase, for example:*
 - Too high load capacity
 - Liquid in the compressor
 - Slugging in pipelines
 - Loss of lubrication due to the dilution of oil

- *Hazards escape of refrigerants such as:*
 - Explosion
 - Fire
 - Toxic
 - Frostbite of Skin
 - Panic
 - Thinning of the ozone layer
 - Greenhouse effect

- *Hazards by moving parts of machinery such as:*
 - Injury
 - Hearing damage due to excessive noise
 - Damage by vibration

- Hazards which are valid for all installations including:
 - Excessive temperatures
 - Incorrect operation
 - Reduction in mechanical strength due to corrosion

- Reduction in mechanical strength due to erosion
- Reduction in mechanical strength by thermal claimed stress
- Reduction in mechanical strength due to vibrations

2.9.5 Noise

The A-assessed equivalent extended sound pressure level at the operating workplace of the KÜHLMOBIL during normal operation at a distance of 2 m without reflection.

There may be a high sound pressure level in the vicinity of the KÜHLMOBIL depending on local conditions.

This increased sound pressure level can be caused by neighbouring devices/machines, and can cause hardness of hearing.

3 Product Description



Note!

Some of the equipment described in these operating instructions is only available on certain the KÜHLMOBIL models.

*This optional equipment is identified by an *.*

3.1 Explanation of Terminology

Secondary:

Water supply: Cold water leaving the KÜHLMOBIL or entering the device to be cooled.


Water return: Heated water leaving the device to be cooled or entering the KÜHLMOBIL.

primary:

Water inlet: Cold water leaving the KÜHLMOBIL or entering the device to be cooled.

Water outlet: Heated water leaving the device to be cooled or entering the KÜHLMOBIL.

3.2 Overview

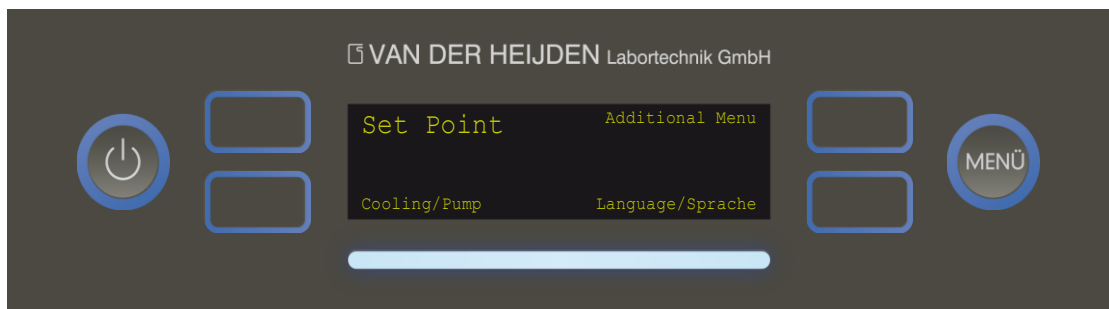

Operating elements on the touch-pad control panel



Figure main switch on the front panel door (at large Chillers)

3.3 Functional Instructions

1. The appliance is customized for the customer's voltages. Therefore the feeder had to be connected before start up. A clamp box for this is mounted on the back side of the chiller. The connection has only done by an expert!



When connecting the supply voltage variation necessarily the note to damage the cooling cell to avoid.

1. When connecting the supply voltage to the three-phase system, the phase sequence must be observed. The KÜHLMOBIL is wired for a clockwise rotating field and must be connected accordingly.

The power supply for a KÜHLMOBIL is set up so that it can be, abge-switched independently of the power supply for other electrical equipment, in particular for lighting-processing equipment, ventilation systems, alarm and other safety equipment

The appliance's coolant circuit is designed as an „open“ system. The feed pump draws the cold coolant out of the container and pumps via the coolant supply through the unit to be cooled.

The heated coolant flows back through the coolant return and into the coolant container, where it is cooled and begins the circuit again.
2. If faults The coolant circuit of the device is designed as an "open" system The feed pump sucks the cooling liquid from the container, promoting it across the water flow through the unit to be cooled.
3. The high pressure switch protects the compressor in the refrigeration cycle devices (type LG / WG) from overpressure If it is no longer guaranteed that the superheated refrigerant in the condenser can be cooled sufficiently, triggers the pressure switch before dangerous damage may occur.

The fault is indicated on the display and the indicator light turns red In addition, an audible alarm sounds.

Reasons may be: eg heavy contamination of the cooling fins of the condenser or excessive cooling liquid return temperature.
4. The low pressure switch protects the compressor in the refrigeration cycle devices (type LG / WG) before vacuum If, for eg, leakage caused by material fatigue in the refrigeration system and escaping refrigerant, the low pressure switch switches off the KÜHLMOBIL before the system is drained, the air could enter

In addition, the low pressure switch is used as anti-freeze for the plate heat exchanger (evaporator)

The error message is displayed, the indicator light flashes red and an audible signal
5. If errors occur in the refrigerant circuit, the feed pump switches off in Kühlflüssigkeitskreislauf The corresponding fault appears on the display and the indicator light turns red In addition, an audible alarm sounds.
6. Once the coolant in the tank falls below the minimum level, coolant lack appears on the display, the status light flashes amber and an audible warning signal

The KÜHLMOBIL switches off in this state

Only coolant to be added to the coolant reservoir before continuing operation of KÜHLMOBIL s needs
7. The actual temperature of the cooling liquid in the cooling circuit is controlled by the electrical control At the film on the front panel of KÜHLMOBIL s the actual and the desired temperature is displayed

The actual value is the coolant exit temperature on the To-circuit designated "water flow" on KÜHLMOBIL.

8. The coolant outlet temperature is continuously monitored and compared with the set temperature set-point. Once the smallest deviation occurs, the temperature controller attempts to correct this

The correction value is (type LG / WG) immediately passed with compressor devices to a corresponding regulating proportional valve in the refrigerant circuit to compensate for the disturbance

The precise control mode can also be achieved by means of tank heater Through permanent control of heating, a high temperature accuracy can be achieved Since the scheme with heater temperature accuracy (0.1K / min) is required, the cooling does not switch off, but runs continuously

The correction value is immediately passed on in water-water units at a correspondingly-regulate the multiway cooling water control valve on the primary side to compensate for the disturbance.

9. Also for devices without power control the coolant outlet temperature is continuously monitored and compared with the set temperature setpoint The scheme utilizes a symmetric hysteresis In this case, the cooling is switched off at the lower limit and switched on again at the upper limit

The fault is displayed and the status indicator flashes red. In addition, an audible warning signal sounds.

10. Because differently sized heat loads can occur, this load change must be continuously monitored and balanced The need exists because high temperature constancy is required.

In large heat load changes is by the power control a correction time of approximately 10 - requires 12 minutes. Thereafter, the corresponding temperature of 0.1 K stability is reached again.

Attention⇒! If too low a heat load supplied, the constancy of the flow temperature is no longer guaranteed by the power control

The result is a reduction in the supply temperature down to the cut-off point of the cooling. However, the pump continues to operate

11. The target temperature is set at the factory This temperature is sufficient for most applications.

In any case, an adjustment of the desired temperature within the specified range is possible Should it happen that in spite of the set coolant flow temperature a tie-Fere actual temperature is displayed this is a result of too low thermal load (to be cooled equipment is too little heat to the coolant from)

In this case, it is not possible to compensate for the full power range. On the consistency of the actual temperature in the water supply, however, this phenomenon has no effect In order to achieve a high temperature stability, it must be ensured that the KÜHLMOBIL operates continuously

The continuous operation is achieved by continuous power adjustment In order not to unnecessarily consume energy, was deliberately omitted on a counter heating

The emitted outward cooling capacity is thus reduced.

12. If the specification is too far below switches the low-temperature protection (not the low temperature alarm) the cooling off The circulation pump continues to run and it appears no alarm

Once the desired temperature in the refrigeration cycle is reached again, the refrigerant circuit is switches on again.

13. If the coolant temperature flow temperature to set temperature alarm turns the unit off completely

The corresponding fault (high temperature alarm or low temperature alarm) appears on the display In addition, an audible alarm sounds and the indicator light flashes red.

3.4 Accessories and Options

3.4.1 Floating Contact

The floating contact is closed when the KÜHLMOBIL is in normal operation. The floating contact opens if the KÜHLMOBIL malfunctions / breaks down.

Especially if the device to be cooled is highly temperature-sensitive, it can be shut down immediately or the shut-down procedure can be initiated in order to prevent further damage.

The floating contact is connected to the KÜHLMOBIL via a two-pin plug-in contact and this work should be done by a trained expert.

The plug is supplied with the KÜHLMOBIL.

The rated voltage for the floating contact is 250V AC / DC. The rated current is 10A AC / DC.

These data are maximum values and are not to be exceeded. (see technical data)

3.4.2 Remote Control

In the case of a KÜHLMOBIL with remote control, the touch-pad control panel is not mounted on the device itself, but rather in separate console housing with a correspondingly long connection cable to the KÜHLMOBIL.

The same functions can be performed on the remote control as can otherwise be performed directly on the KÜHLMOBIL.

A total cable length of up to 40 metres is possible.

3.4.3 Operating Hours Counter

The operating hours counter indicates the pure operating time excluding standby time.

The operating hours counter is mounted beside the touch-pad control panel on the front of the KÜHLMOBIL.

3.4.4 Flow indicator

With the flow indicator the current water flow in the cooling water return is shown. The required flow rate should be by the rear bypass valve.

With several flow indicators the flow is shown separately for each circuit.

Since the flow is available in different versions, is here no application.

3.4.5 Flow control unit

The current coolant flow in the secondary water return is measured with the flow control unit and displayed on the front display.

The required flow rate is adjustable through the bypass valve

The flow meter is internally installed in the device and to adjust only by a professional using.

Note ⇒ the currently funded amount of coolant is visible on the front panel display (temperature).

The display shows the amount in L/min.

3.4.6 External switch on

With the external involvement of the KÜHLMOBIL can from any other location switched on or off. If the external activation via the rear plug contact is activated, the KÜHLMOBIL starts automatically.

The external terminal supplies a voltage a potential free contact. When connecting the switching voltage must be observed to prevent damage to the unit.

The connection of external involvement in the KÜHLMOBIL is via a two-pole plug-in contact and should be performed by a trained professional, technician.

The connector is supplied with the KÜHLMOBIL s.

After starting up the KÜHLMOBIL external intervention must be connected. Only then can the KÜHLMOBIL via the front panel turn on and operate.

Is the plug contact closed with active cooling, the control panel is disabled and the display shows:

„external switch aktiv“

If the plug-in contact closed before switching on the pump and the cooling, the control panel is disabled and the display shows

„external switch aktiv – “please disconnect”

and the status indicator flashes green.



3.4.7 Flow monitor

The flow monitor in the secondary water supply switched on the device at low flow from through the device being cooled.

The flow monitor in the secondary water supply shut down the unit in case of insufficient the pump. It serves to protect the pump.

The flow monitor in the secondary return water behind the internal plate heat exchanger shear turns the unit at low flow from through the heat exchanger. He serves as an ice protection.

The shutdown value set at the factory. Subsequent adjustment is possible at any time. The flow switch is installed in the device and to adjust only by an expert.

Since the flow monitor is available in several versions, is omitted here an illustration.

3.4.8 Time Relay for Run-On

The time relay is not visibly installed inside the housing. It ensures that the circulating pump runs on for a preset period of time (e.g. to remove the residual heat from the device to be cooled) after the KÜHLMOBIL is switched to standby.

The run-on time (specified by the customer) is usually set prior to delivery to the customer and is subsequently permanently set.

3.4.9 Overflow

For security reasons KÜHLMOBILs are additionally equipped with automatic refilling with an overflow. This serves as a protection to prevent when a defective solenoid valve Automatic refill overflowing of water tanks.

The terminal "Overflow" is back out on the heat out mobile. This is a hose.

Attention: ➡ The overflow has to be connected and be guided with gradient into an unpressurized flow!

The overflow not constrict under any circumstances!

3.4.10 Automatically re-fill

An automatic refill within the KÜHLMOBIL ensures that the water tank of the secondary circuit is always filled optimal.

If a water loss can occur, a level switch is activated, which in turn activated a magnetic valve allowing that water flow into the water tank.

A prerequisite for this function is, however, a constant supply of connection "refill" water.

The automatic refill with water can cause the flow at temperatures below +10°C is not sufficient antifreeze.



In the event of a leak in the cooling system can lead to an uncontrolled automatic refill of the water outlet.

For this reason, the corresponding area on the KÜHLMOBIL, at the hose connections and at the device which have to be cooled must be monitored with moisture sensors.

By a leakage to avoid further damages, it must be ensured that all water supplies and the KÜHLMOBIL shut off immediately and be protected against restart!

3.4.11 Bypass

For version of the cooling circuit with the bypass regulating valve required for setting is available on the rear of the chassis. This is the possibility of regulating the cooling fluid flow on the flow rate and pressure so as to influence the flow rate of the circulation pump.

It should be taken to ensure that on opening of the bypass valve of the cooling liquid pressure is sufficient in any case.

The supply pressure is back on the pressure gauge (option) to read and can be adjusted to the desired value. Clockwise that it is not means pressure increase and counter-clockwise pressure reduction.

Please note unnecessarily high pressures are set.

The lower the pressure is displayed on the pressure gauge, with adequate number of liters, the better for the coolant pump.

When **automatic bypass** the regulation of the flow rate is automatic. From the main water flow not currently required coolant quantity is diverted from the water supply and fed directly through the bypass into the return water from the main flow.

The automatic bypass valve is factory-set to an opening pressure. An adjustment is possible at any time.

3.4.12 Manometer

The manometer indicates the coolant pressure in the cooling water return.

In the case of the KÜHLMOBIL with a bypass, the supply pressure is indicated on the pressure gauge.

No picture is supplied here, as the pressure gauge is available in different versions.

3.4.13 Dirt Filter

The dirt filter is built into the cooling water supply. It removes any dirt particles from the cooling water before it enters the device to be cooled.

The fine filter is installed in the secondary coolant supply. It cleans the cooling water before it enters the unit which is to be cooled by any existing dirt particles.

No picture is supplied here, as the dirt filter is available in different versions.

3.4.14 Fine Filter

The fine filter is installed in the primary water inlet or the secondary water supply. It cleans the cooling water before entering the KÜHLMOBIL or before entering the cow to Loin device from any existing debris.

Since the fine filter is available in several versions, is omitted here an illustration.

3.4.15 Clean dirt filter

With decreasing power and / or low flow the strainer must be cleaned in the water inlet (for water-cooled units) or water flow.

To clean / replace the strainer, the device must shut down, drained the coolant and the primary coolant are blocked in water cooled KÜHLMOBIL.

Only after this the strainer must be opened.

Can the dirt-filter be washed out he must not be renewed - cleaning ranges.

3.4.16 Clean fine filter

With decreasing power and / or low flow of the fine filter must be cleaned (for water-cooled units) and the water flow in the water inlet.

To clean / replace the filter, the device must shut down, drained the coolant and the primary coolant are blocked in water cooled KÜHLMOBIL. Only after that it's allowed to open the filter.

Can the fine-filter be washed out he must not be renewed - cleaning is enough.

3.4.17 Interface RS232

The interface adapter is used for remote control and VDH_SER1 error reading via an RS232 interface of KÜHLMOBIL that are equipped with the regulator VDHR201. To avoid interference, the interface of the controller is electrically isolated.

Communication with the controller takes place with a simple textorientierten protocol. It can be integrated into existing programs.

To test the functionality, when using a PC with Windows - the operating system, for example the accessory HyperTerminal program may be used.

Interface parameters: 9600 baud, 8 data bits, 1 stop bit, no parity, no flow control.

3.4.18 Geodetic height difference

The geodetic height difference is an option which in the coolant supply a non-return valve and a solenoid valve is built into the coolant return. The option must be built if the KÜHLMOBIL lower the cooling unit is at least 3 m.

The valve prevents that collapses after shutdown of the pump(s) of the KÜHLMOBIL the water and overflows coolant tank.

3.4.19 Tank drain

When connecting the tank emptying is a ball valve with hose.

The tank drain should be connected with a slope at an unpressurized drain or another vessel.

Notes!

The connection of the tank emptying are back on KÜHLMOBIL and is angled downward.

- Ensure that KÜHLMOBIL set up so that the water hose to be connected between KÜHLMOBIL and outflow can not be exceeded.
- Ensure that KÜHLMOBIL set up so that the water hose to be connected between the drain and KÜHLMOBIL not located in passageways.

3.4.20 Automatic Switchover to Water Pipe in the Event of a Breakdown

In the event of a breakdown of the KÜHLMOBIL or a power failure the cooling-water supply and return are switched over to the on-site water supply pipe.

This changeover remains active as long as the malfunction or power failure remains unremedied.

Two connection nozzles protrude from the back of the KÜHLMOBIL to which the pipes for water supply and drainage can be connected.



In the event of a leak in the cooling system can lead to an uncontrolled automatic refill of the water outlet.

For this reason, the corresponding area on the KÜHLMOBIL, at the hose connections and at the device which has to be cooled must be monitored with moisture sensors.

By a leakage to avoid further damages, it must be ensured that all water supplies and the KÜHLMOBIL shut off immediately and be protected against restart!

3.5 Protective Equipment

The KÜHLMOBIL is fitted with protective equipment so that there are no hazards to the safety and health of the operator or of third parties when the unit is being used correctly.

The KÜHLMOBIL is designed and built in accordance with state-of-the-art technology and recognised safety rules.

In order to guarantee safe operation for personnel and the environment, the precautions described below have been introduced and safety equipment has been installed.

Unauthorised removal or bypassing of protective equipment constitutes an indictable action.

In the event of damage/injury, any liability claims will be forfeited.

All components, which can become hot or cold, are accommodated inside the housing, and are also secured by means of additional protective panels/grids, which can only be removed using tools.

In the event of a power failure or a reported malfunction, all elements of the control system revert to a safe status for operator, KÜHLMOBIL and the environment.

All voltage-carrying components are secured (insulated) against contact and installed at sufficient spark-over distances.

Voltage-carrying modules are only installed inside the housing, and can only be opened using tools.

All electrical components bear the CE identification for low voltage and/or EMC.

To earth the KÜHLMOBIL, potential equalization ($\varnothing > 1.5 \text{ mm}^2$) is installed for all conductive components.

The KÜHLMOBIL is designed in accordance with protection class IP 52.

All control-unit connections are clearly marked.

Both the electrical and hydraulic circuit diagrams can be found in the appendix together with their accompanying spare-parts lists.

4 Transportation, Storage and Commissioning

4.1 Transportation

The KÜHLMOBIL is delivered in a vertical position in a wooden crate. The following options are available for transporting the wooden crate :

- Forklift and
- pallet truck.



Danger!

Danger to life due to the wooden crate and the KÜHLMOBIL contained inside falling!

As a result of incorrect transportation of the wooden crate the wooden crate and the KÜHLMOBIL contained inside can fall and cause very serious injuries or death.

Use only suitable lifting equipment and securing elements!

The driver must be entitled to drive the forklift/pallet truck.

Follow the following instructions for transporting of the KÜHLMOBIL:

The forklift/pallet truck must be certified for the total weight of the wooden crate plus the KÜHLMOBIL (see Section 9.1, *Technical Specifications* for weight specifications).

The forks of the forklift/pallet truck must be of sufficient length.

Before the KÜHLMOBIL is lifted, all persons must vacate the working area of the forklift/pallet truck.

The KÜHLMOBIL is a sensitive device and must be transported with corresponding caution in its wooden crate as far as the erection site.

Avoid the wooden crate hitting the ground abruptly when setting it down (e.g. on the loading surface of a transportation vehicle).

The wooden crate is to be secured against sliding and slipping down in the transportation vehicle.

The wooden crate is to be secured against tipping over during transportation.

The wooden crate absolutely must be transported and stored in dry conditions.

Arrows printed on the transportation crates must point upwards.

Further instructions on the packaging must also be followed!

Adhere to applicable accident-prevention and work-safety regulations.

4.2 Unpacking/Scope of Delivery



Instructions!

Unpack the KÜHLMOBIL immediately after delivery.

Upon delivery, be sure to look out for sever damage or damage due to incorrect transportation!

Proceed as follows when unpacking the KÜHLMOBIL:

- Slacken all bolts on the lid and on the front of the wooden crate
- Remove the cover and front panel.
- Loosen the screws with those the internal framework for the adjustment of KÜHLMOBIL are screwed on.
- Pull to the frameworks upward from the transportation crate.
- Remove the upper Styrofoam panels from the KÜHLMOBIL (if present).
- Remove the KÜHLMOBIL using appropriate equipment (forklift/pallet truck) from the wooden crate ⇒ See 4.1.
- Remove the Styrofoam elements from the underside of the KÜHLMOBIL (if present) as soon as they are accessible.
- Check that the delivery is complete on the basis of the *delivery docket*.
- Check the entire delivery for external damage that might, for example, have been caused during transportation.
- If transportation damage is discovered, contact the supplier/delivery company immediately.
- **Complaints at a later date will not be accepted!**
- If parts are missing or have been delivered incorrectly, contact. **Van der Heijden Labortechnik GmbH**.

4.3 Put Up



Note!

Unpack the KÜHLMOBIL immediately after delivery.

When erecting the unit, ensure that there is sufficient room for manoeuvre while working.

Do not erect the KÜHLMOBIL with the air-intake side directly opposite heating elements/heaters.

Proceed as follows when put up the KÜHLMOBIL:

- Roll the KÜHLMOBIL carefully and cautiously to the erection site.
- Erect the KÜHLMOBIL on an even, solid base.
- Secure the KÜHLMOBIL against rolling away unintentionally by pushing down the locking brake on the rollers.



Instructions!

The air flowing backwards through the condenser must not be hindered due to an insufficient gap. The hot air must be able to flow away upwards unhindered.

When erecting the unit, it must be remembered that the KÜHLMOBIL releases its cooling output into the surrounding room as heat.

The produced heat is increased even more by the drive output of the cooling aggregate.

The resulting increased room temperature would reduce the cooling capacity of the KÜHLMOBIL.

For this reason sufficient air conditioning of the room in which the unit is erected must be guaranteed.

- Position the KÜHLMOBIL so that the wall is at least 80 cm away from the front of the unit to absorb the cool air.
- Position the KÜHLMOBIL so that the wall is at least 80 cm away from the back of the unit.
- Please take care that the heated air could not be absorbed on the front side.
- When the tank is filled with water please do not move the KÜHLMOBIL anymore.

4.4 Storage

The KÜHLMOBIL must be stored in a vibration-free, dry and as dust-free a location as possible. It must not be stored outside locked rooms.

The air temperature must be within a range of + 5°C and + 32°C.

Relative humidity must not exceed 85 %. Condensation of air-borne water vapour on the surfaces of the KÜHLMOBIL absolutely must be avoided.

4.5 Starting Up the KÜHLMOBIL



Danger of death due to electrocution!

There is a risk of electrocution if incorrect work is carried out on electrical components!

Work on electrical systems may only be performed by specialist electricians.

There is a danger of electrocution if improper work is carried out on electrical components!

For devices with voltage changeover applies:

Before the first terminal of KÜHLMOBIL to the supply voltage of the voltage selector switch must be set by authorized technicians for installation.

Please note that the factory setting! (see specifications).

Subsequent adjustment at a later time, for example, by the operator shall be prohibited.

Only in this way damage to the KÜHLMOBIL can be avoided.



Caution!

Irreparable damage to the cooling circuit is possible!

The oil in the coolant circuit of the KÜHLMOBIL must first accumulate before switching on.

Do not switch on the KÜHLMOBIL for at least one hour after the unit has been positioned at the erection site.

4.5.1 Bleeding the Circulating Pump

The circulating pump (feed pump) of the coolant circuit must be bled in order to enable the coolant to be pumped:

- Carefully fold the cover of KÜHLMOBIL backwards. (Up to size 08)
- Unscrew the stopper / cap from the coolant tank.
- Fill the underlying coolant reservoir to below the edge of the pool with coolant. (See max. Level specified in the Appendix).
- Connect a piece of hose to the water supply.

4.5.2 Cooling-Water Hoses



Instructions!

Both cooling liquid connections are quick-release fastener connections with plug-in nozzles.

The supply pressure is to be taken into account when choosing the cooling-water hoses for *water supply and return*.

Terminology definitions:

Water supply Cold water leaving the KÜHLMOBIL or entering the devices to be cooled.

Water return: Heated water entering the KÜHLMOBIL or leaving the devices to be cooled.

- Attach the cooling water hoses for water supply and return between the device to be cooled and the KÜHLMOBIL.
- First connect the cooling water hoses for *water supply and return* to the devices to be cooled (instructions in the operating instructions for the respective device).



Instructions!

The hose nozzles are located securely mounted on the back side of the KÜHLMOBIL.

- Ensure that the KÜHLMOBIL is erected such that the water hoses to be connected between the KÜHLMOBIL and the devices to be cooled do not lie in walkways.
- Ensure that the KÜHLMOBIL is erected such that the water hoses to be connected between the KÜHLMOBIL and the devices to be cooled do not be stood on.



Note!

Temperatures in the cooling-water return (e.g. due to connection of additional devices to be cooled) in excess of + 30°C are to be prevented in order to prevent overloading of or damage to the compressor.

- Secure the cooling water hoses to the respective connection at the rear of the KÜHLMOBIL.
- Do not confuse the *water supply* with the *water return*. If in doubt it is better to double check

4.5.3 Primary

When choosing the coolant hoses for water inlet and outlet of the pressure range is observed.
(see specifications)

To establish the terminology:

Water inlet: Outlet of cold water from the “in-house circulation” or entry into the KÜHLMOBIL.

Water outlet: Outlet of the heated water from the KÜHLMOBIL or entry into the “in-house circulation”.

- Route the coolant hoses for water inlet and outlet between the in-house circuit and the KÜHLMOBIL.
- Now close the coolant hoses for water inlet and outlet of the primary circuit of the KÜHLMOBIL.

NOTE!

The connectors are located at the back on KÜHLMOBIL and are angled downward.

- Make sure that the KÜHLMOBIL install so, that the water hoses to be connected between KÜHLMOBIL and the in-house circulation are not across passageways
- Make sure that KÜHLMOBIL install so that the water hoses to be connected between KÜHLMOBIL and the in-house circulation can not be exceeded.



In the event of a leak in the cooling circuit can occur through the primary cooling water to an uncontrolled water outlet.

For this reason, must the corresponding area on KÜHLMOBIL, in the manifold and on the device to be cooled, monitored with moisture sensors.

In order for a leak to avoid further damage must be ensured that all water supplies and the KÜHLMOBIL immediately switched off and on again-turn, be secured.

4.5.4 Floating Contact

The floating contact is closed during normal operation of the KÜHLMOBIL.

The floating contact opens in the event of a malfunction/breakdown of the KÜHLMOBIL.

In particular if the device to be cooled is highly temperature sensitive, it can be shut down immediately by the switching contact or the shut-down procedure can be initiated in order to prevent further damage.

The floating contact is connected to the KÜHLMOBIL by means of a two-pin plug-in contact.

The plug is supplied with the KÜHLMOBIL.

The floating contact is connected to the back of the KÜHLMOBIL.

- Have the two-core wire for the floating contact attached and subsequently tapped onto the plug by an electrical expert.
- It is essential that the plug then be connected and secured.



Note!

If the KÜHLMOBIL is supplied with both a floating contact and an external switch-on unit * (see below), then the two individual connections are brought together in a five-pin plug-in contact, or equipped with two plug-in contacts with different (unmistakable) plug-in connections!

4.5.5 External Switch-On Control

With the external switch-on control the KÜHLMOBIL can be switched on and off from any other desired location – e.g. from the electron microscope.

The external switch-on control is connected to the KÜHLMOBIL by means of a two-pin plug-in contact. The plug is included with the KÜHLMOBIL.

The external switch-on control is connected at the rear of the KÜHLMOBIL.

- Have the cable for the external switch-on control attached and subsequently tapped onto the plug by an electrical expert.



Note!

If the KÜHLMOBIL is to be connected to an external switch-on control, the contact (plug-in contact on the KÜHLMOBIL) to the KÜHLMOBIL must be due to the configuration connected. Before activation of the external switch on control the KÜHLMOBIL must be start up (pump and cooling).

If the external switch on control is activated before start up the display shows:

External switch on control active – please cut



The status lights are green – flashing.

- If the specification on the plug is *floating free open* ☞ „external switch on control active“ the external switch on control must be connected **before** start up the KÜHLMOBIL.
- If the specification on the plug *free from tension* ☞ „external switch on control active“ the external switch on control must be connected **before** start up the KÜHLMOBIL.
- If the specification on the plug *floating free closed* ☞ „external switch on control active“ the external switch on control must be connected **after** start up the KÜHLMOBIL.
- If the specification on the plug *voltage on the plug* ☞ „external switch on control active“ the external switch on control must be connected **after** start up the KÜHLMOBIL.



Note!

If the KÜHLMOBIL is supplied with both an external switch-on unit (see above) and a floating contact, then the two individual connections are brought together in a five-pin plug-in contact, or equipped with two plug-in contacts with different (unmistakable) plug-in connections!

4.5.6 Remote Control

In the case of a KÜHLMOBIL with remote control, the touch-pad control panel is not mounted on the device itself, but rather in separate console housing with an appropriately long connector cable to the KÜHLMOBIL.

The same functions can be performed on the remote control as can otherwise be performed directly at the KÜHLMOBIL.

The remote control is connected to the KÜHLMOBIL via a 4-core shielded cable.

- Have the remote control cable attached by an electrical expert.

5 Touch-Pad Control Panel

The functions of the touch-pad control panel are described in greater detail, and the actual operation of the KÜHLMOBIL is described in the subsequent chapter.

The KÜHLMOBIL is operated from the touch-pad control panel.

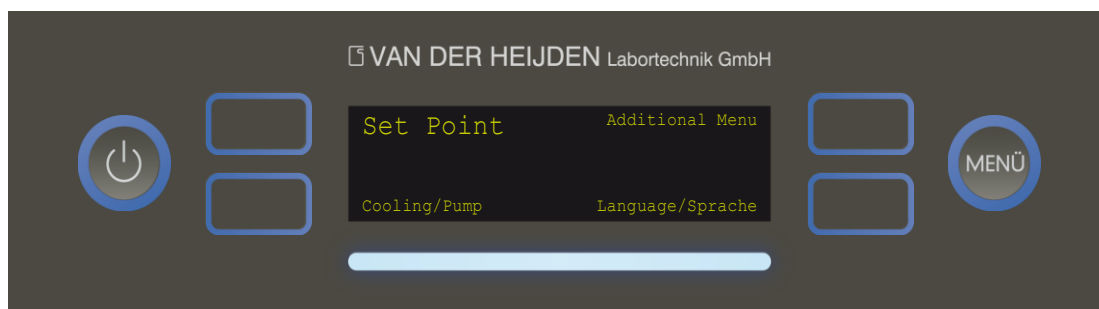
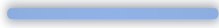
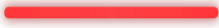
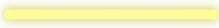


Fig. 1: *Touch-pad control panel*

Operating elements on the touch-pad control panel of the KÜHLMOBIL:

5.1 Operating Features

Illustration	Control Functions
	On/Off Button: the KÜHLMOBIL (Mobile Chiller) is switched 'ON' and 'OFF' via this control. Once the supply voltage is connected to the control facility, this button is illuminated 'blue'. <i>Exception:</i> the button is only illuminated 'blue' when the power supply is activated externally, but has no other function than a switch.
	The 'Menü' Button enables a changeover between the Operating Display, the Switch-on- and Parameterage Display and the Information Display. • The <u>Operating Display</u> shows all the operating information - such as the operating state and temperatures - necessary for a disruption-free functioning of the appliance. • The <u>Switch-on- and Parameterage Display</u> enables settings to be made for the nominal parameters and the chiller unit. • The <u>Information Display</u> enables call-ups to be accessed of parameterage, the fault memory and the appliance information.
	The various Menu levels enable functions to be displayed and/or set by up to four function buttons.
	Display state blue: - the appliance is connected to the power supply. - refrigeration and pumps are running. - no disruption message is shown.
	Display state green: - flashing green illumination when the regulator unit is on stand-by - o only the control facility is activated in stand-by. - o all menu functions can be accessed. - the display state switches to continuous green illumination once a pump is running.
	Display state red: - disruption message shown – an acoustic tone is also emitted. - the disruption cause can be shown on the display. - the acoustic tone is switched-off by calling-up the disruption cause via the function button. - the display repeats flashing red. - follow the menu items to rectify the disruption.
	Display state yellow: - the display flashes yellow when the fluid level is low. - o proceed as with display state red. - once the disruption (display state illumination flashes red or yellow) is rectified, the display state switches to continuous yellow illumination. - o the rectification of a disruption should be acknowledged via the function button. - o the appliance switches back automatically to the general operational state as before the disruption (display state turns to blue illumination).

6 Operating the Appliance

6.1 General


6.1.1 Operating Conditions

The KÜHLMOBIL (Mobile Chiller) appliance is best operated in a dry and dust-free location with no exposure to vibration. The KÜHLMOBIL (Mobile Chiller) should not be operated outside of enclosed roomspace (except when the appliance is fitted for outside use). The ambient air temperature should be in a range of +5°C to +32°C. The relative ambient humidity should not exceed 85%. Any condensation from ambient humidity on the casing surfaces of the KÜHLMOBIL (Mobile Chiller) should absolutely be avoided.

6.2 Switching 'ON'


The KÜHLMOBIL (Mobile Chiller) appliance is brought to the stand-by state as follows:

- Connect the KÜHLMOBIL (Mobile Chiller) with the power supply cable provided via the rear refrigeration unit connector to the appropriate power supply.

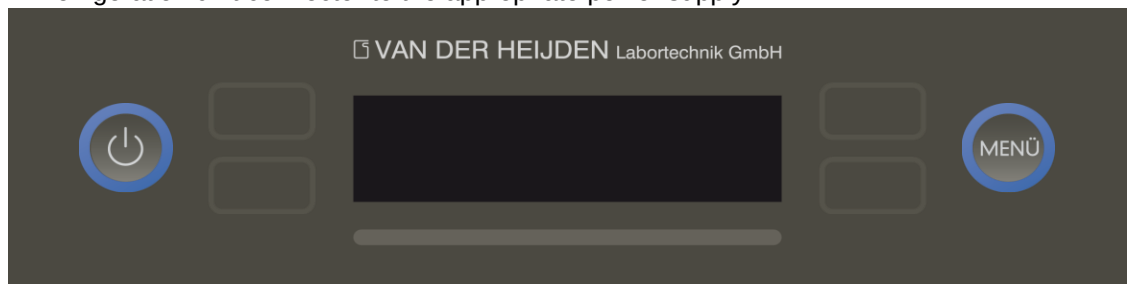

-Fig. 2: the KÜHLMOBIL (Mobile Chiller) appliance when 'OFF'.

- The operating button and the menu button are illuminated blue and switched to functional.
- Depress the  button. The nominal parameter and the operating status of pump and refrigeration are then displayed.

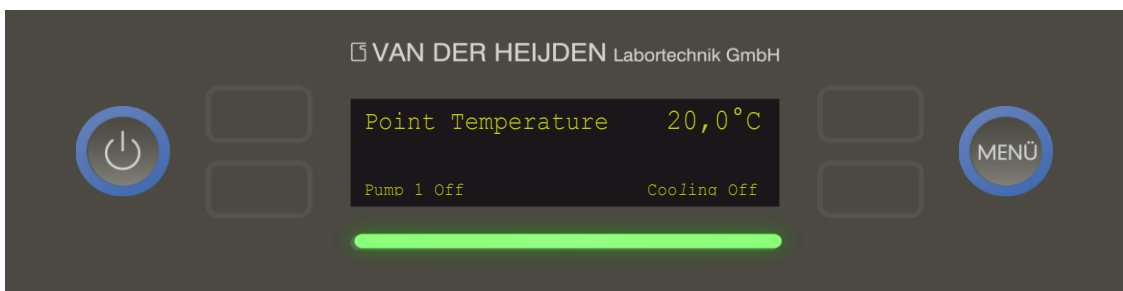

-Fig. 3: 'stand-by'.

- The display status below the display flashes green.




N.B.

Should a disruption message be displayed, or an acoustic warning tone be emitted, refer to the instructions in Section 6.8. The display status flashes.

- Try opening the faucets, which you have installed in the chiller fluid lead between the KÜHLMOBIL (Mobile Chiller) appliance and the refrigerating unit.

6.3 Parametering and Switching 'ON'

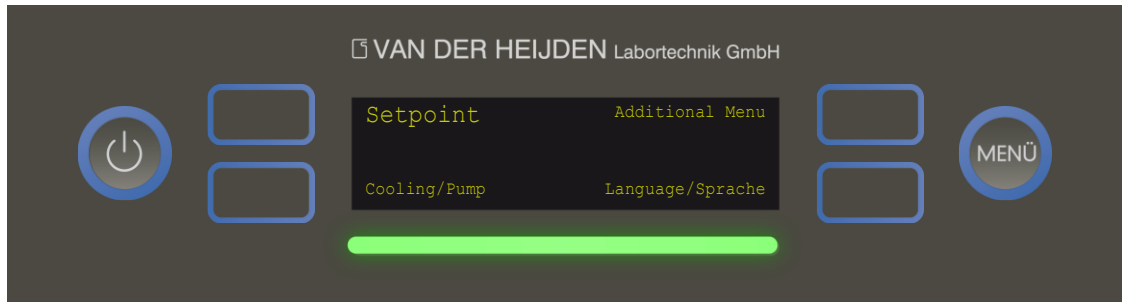


-Fig. 4: 'Parametering and switching 'ON''-


- Depress the menu  button and access the parametering- and switching-on display (see Fig. 4)
- Set the desired language with the  function button.
- By depressing the function button  set the required nominal temperature.
Setting range: The temperature setpoint can be adjusted within limits. (see specifications)
- The KÜHLMOBIL (Mobile Chiller) is now set to the desired user mode and can now be switched-on via the function button: 'pump/chilling'.

6.4 Switching 'ON'

- Depress the function button 'pump/chilling'.

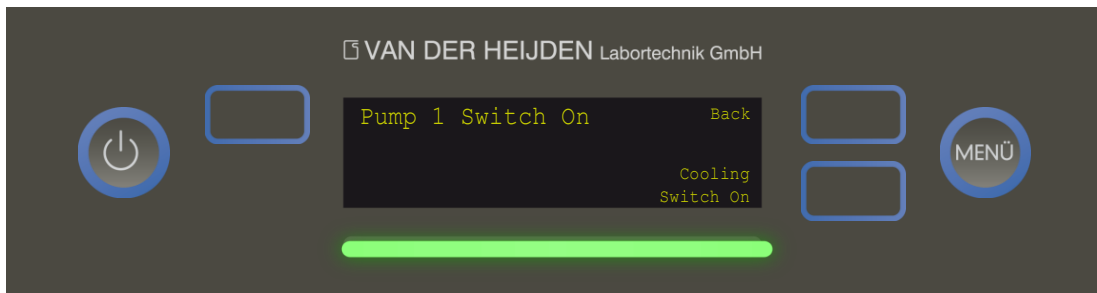

- Depress the function button to switch-on the pump.
- Then, re-depress the function button to switch-on the chilling facility.
- The display state is illuminated blue.

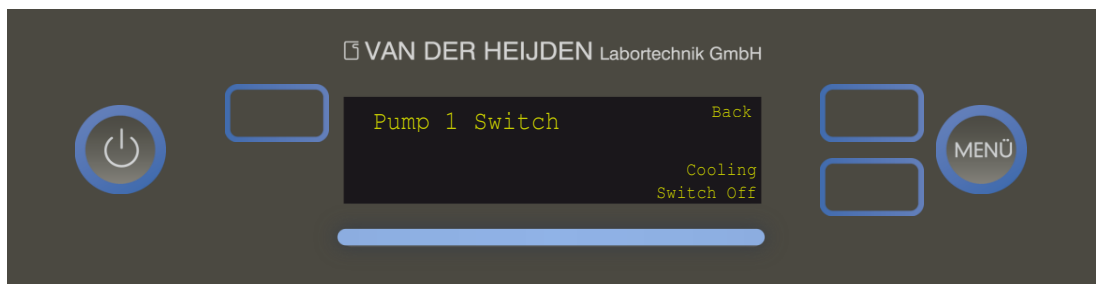

- Use the function button 'return' to re-access the menu: 'parametering and switching 'ON''.

- Use the function button  or wait a dwell time of 60 seconds, when the general operating state is displayed with all relevant information.
- In addition to the references to pump and chilling, a symbol is also activated (illustrated here by a circle), which describes the function.
- Only now connect the plug to the socket of the external switching-on facility (this is an option) on the rear of the KÜHLMOBIL (Mobile Chiller) appliance. When the connection is correctly made, then the display will show the message: 'external switching-on facility'.

i

N.B.

The operating buttons at the KÜHLMOBIL (Mobile Chiller) appliance will be deactivated by connecting up the external switching-on facility.

This will prevent a third party from being able to access the operating panel directly at the KÜHLMOBIL (Mobile Chiller) appliance.

i

N.B.

Only the above described switching-on sequence will ensure, that the aforementioned settings entered for operating the chiller will be saved in the control facility.



-Fig. 5: when in the 'general operating state'

The KÜHLMOBIL (Mobile Chiller) is now set up for the continuous operating mode.

6.5 Continuous Operating Mode

The following routines are to be carried out regularly during the running of the continuous operating mode:

- Regularly check the fluid level in the chiller fluid tank of the KÜHLMOBIL (Mobile Chiller) appliance.
- In case of any obviously loss of chiller fluid within a short period of time, then the circulation system between the KÜHLMOBIL (Mobile Chiller) and the appliance to be chilled, should be checked and any leakage rectified.
- Regularly check and clean the air inlet of the condenser with oil free Compressed air or nitrogen.
- Regularly check the state of the chiller fluid (for algae and bacterial slime).
- When a problem with algae possibly add Thermo Clean DC (see accessories and spare parts list)

6.6 Switching 'OFF'


i

N.B.


The KÜHLMOBIL (Mobile Chiller) appliance is designed for continuous running. The appliance can thus be operated over a longer period of time without interruption.

To conserve energy, the appliance should be switched-off during longer time periods of inactivity.

- You should ensure, that when switching-off the KÜHLMOBIL (Mobile Chiller) appliance, the operations are not endangered of the connected-up appliance to be chilled.
- Deactivate the external switching-on facility via the electron microscope, or disconnect the electrical connector plug of the external switching-on facility * on the rear of the KÜHLMOBIL (Mobile Chiller) appliance. The message on the display: 'external switching-on facility' will then extinguish.
- The chiller unit and the pumps are now switched-off. The KÜHLMOBIL (Mobile Chiller) appliance is now switched-over into the 'OFF' mode.

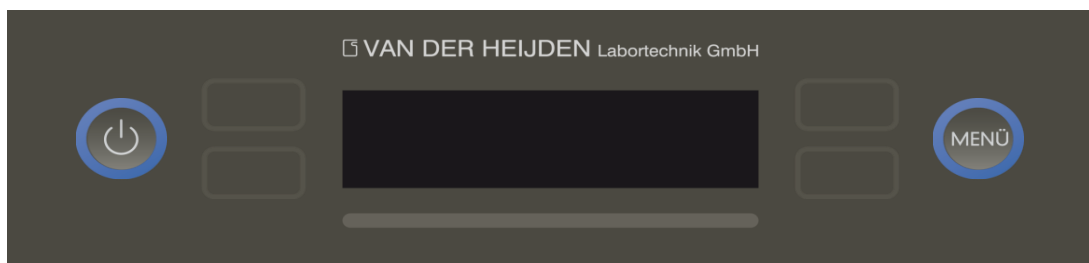

- In order to switch-off the KÜHLMOBIL (Mobile Chiller) without the Option: 'external switching-on facility', depress the button: .
- An additional message will require the procedure to be acknowledged, in order to avoid any inadvertent switching-off of the appliance.
- The chiller and the pumps are now switched-off, and the KÜHLMOBIL (Mobile Chiller) appliance is now switched-over to the 'OFF' mode.

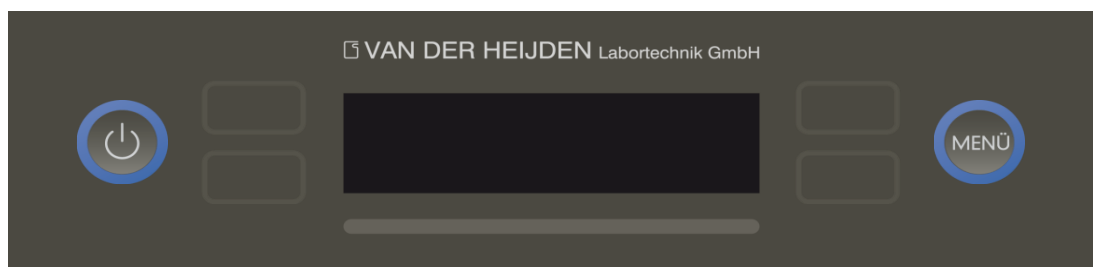

- In order to switch-on the KÜHLMOBIL (Mobile Chiller) again, please refer to Chapter 3.8. Once the button  is depressed, the KÜHLMOBIL (Mobile Chiller) appliance returns to the operating state in which it was set before switching-off.

i


N.B.

In case of an electrical power cut, the following is to be observed:

Once the electrical power supply is restored after a power cut, the KÜHLMOBIL (Mobile Chiller) appliance returns to the operating state, in which it was set before the power cut. No new setting- or control procedures are therefore necessary at the KÜHLMOBIL (Mobile Chiller).

6.7 Information Display

The information display will show parameterage, fault memory, appliance information and display settings when called-up with the function buttons.

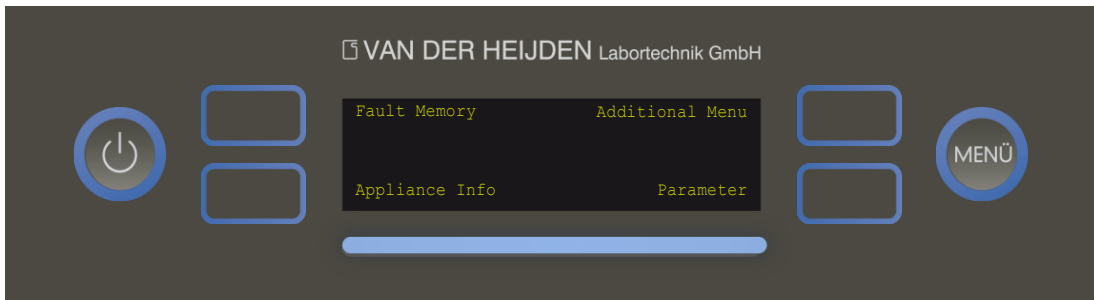

-Fig. 6: Information display-

All control settings can be accessed by the function buttons on the parametering level. Any alterations thereto are not possible in this menu mode.

6.7.1 Parameter

In the parameter level, all settings of the controller can be called up with the function keys. Changes are not possible in this menu


6.7.2 Fault Memory

The fault memory records all disruptions with indication of time-of-day. The memory also records when a fault occurred and when the fault was rectified..

6.7.3 Appliance Data/Details

The date/details of the appliance include the serial number of the KÜHLMOBIL (Mobile Chiller) appliance, among other things, which are the serial number of the control facility, and the postal address and telephone number of the corporation: Heijden Labortechnik GmbH (Inc.).

6.7.4 Display Settings

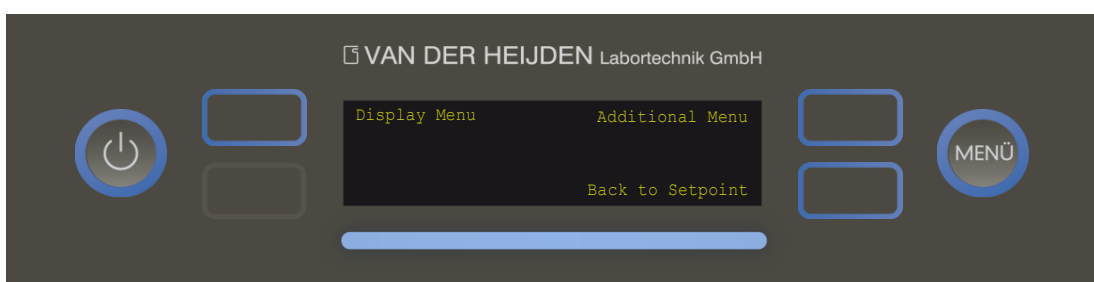

-Fig. 7: Information display (2nd level)

The brightness of the display can be altered by means of the display settings.

6.8 Fault Messages and Disruption Rectification

Fault messages are indicated by the emitting of an acoustic tone, by- a message report on the nominal temperature display, and –a display state red or yellow.

(the relative alarm states are described in the Table down below).



Caution

Irreparable damage can occur to the KÜHLMOBIL (Mobile Chiller) appliance when the cause of a disruption is not immediately rectified.


Please rectify the cause of a disruption with immediate effect.

The acoustic alarm tone is switched-off by depressing the function button.

All disruption messages should be reset by means of the function button after disruption rectification.

6.8.1 Schedule of possible Alarm States:

Disruption: Low fluid level

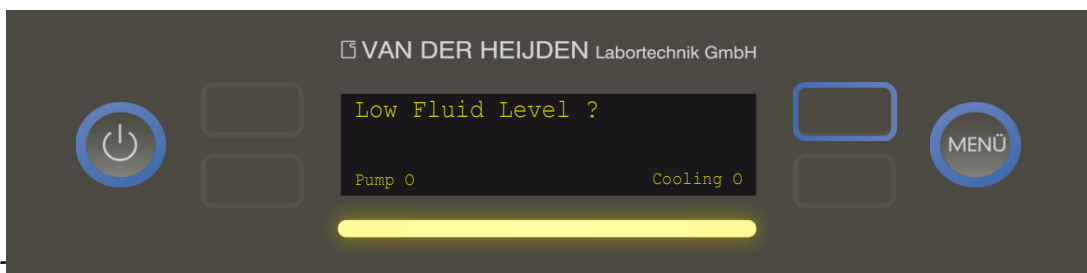

The following changes appear on the display:

- The regulator facility emits a peeping tone and the display status flashes yellow
- The active symbols come to a halt

Disruption rectification:

- The fault should be acknowledged with the function button and the menu accessed for fault rectification.
- Then, follow the instructions.
- Once the fault is rectified, the display status changes to steady illumination and the query of the fault rectification on the display can be acknowledged.
- The display indications return to the disruption-free operational mode.
- The active symbols recommence running and the display status returns to blue.

Other Disruption States:

● <i>Excess temperature alarm</i>
● <i>Deficient temperature alarm</i>
● <i>High pressure disruption – chiller-side</i>
● <i>Through-flow disruption – water-supply-side</i>
● <i>Excess temperature alarm – sensor defect – wiring interruption (sensor)</i>
● <i>Deficient temperature alarm disruption (short circuiting at the sensor)</i>

Other disruptions only occurring individually on technical grounds:

● <i>Low pressure disruption – chiller-side</i>
● <i>Disruption - pump thermal contact</i>
● <i>Alarm – through-flow monitor</i>
● <i>Check fluid level</i>
● <i>Motor contactor relay is triggered</i>
● <i>E4 alarm</i>

i

N.B.

Possible causes and rectification suggestions can be accessed by the function button and found under the menu item 'fault rectification'. Should a fault rectification not be possible, please contact **Van der Heijden Labortechnik GmbH**. Further information can be found in the information display in appliance data/details.

Disruption – Excess temperature rise alarm (example)



-Fig. 9: Disruption – 'Excess temperature rise alarm'

The following changes appear in the display:

- The regulator facility emits a peep tone and the display status flashes red
- The activated symbols come to a halt.

Fault rectification:

- Depress the function button to acknowledge the fault and the menu accessed for fault rectification.
- The, follow the instructions.
- Once the fault is rectified, the display status changes over to continuous yellow and the query of fault rectification can be acknowledged.
- The displays return to the disruption-free operating status.
- The active symbols recommence running and the display status again becomes blue.

6.9 Disruption without Direct Message Reporting

Disruption	Cause	Description / disruption rectification
Pump fails to circulate the refrigerant	<i>Air in the system</i>	– Bleed the circulation pump (see the KÜHLMOBIL (Mobile Chiller) appliance).
	<i>Motor defective</i>	– Contact the manufacturer
	<i>Condenser defective</i>	– Contact the manufacturer
	<i>Long downturn time</i>	– Contact the manufacturer  Caution Damage to the feed pump of the refrigerant circulatory system. The feed pump should be switched-on once weekly for around 5 minutes, with the circulatory system full of refrigerant).
Compressor runs but no chilling effect is produced	<i>Refrigerant level only for units with compressor</i>	– Contact the manufacturer.
	<i>Start relay defective</i>	– Contact the manufacturer.
	<i>Internal capillary cracking</i>	– Contact the manufacturer.
	<i>Capacity regulation without function</i>	(optional) – Contact the manufacturer.

6.10 Display Indications with Appropriate Options

(see specifications)

6.10.1 Two pumps

When an additional chilling circulatory system (optional) is installed with an additional circulation pump (circulation pump '2'), this will be indicated by an additional function button on switching-on.

It is sufficient for one of both the pumps is switched-on (running) for the KÜHLMOBIL (Mobile Chiller) appliance to be functional.

6.10.2 External switching-on

The external switching-on facility is also shown on the display. When this is active, it is not possible to access the KÜHLMOBIL (Mobile Chiller) via the display.

The display is only then re-usable when the external switching-on facility is deactivated.

6.10.3 The run-on time after external switching-on

After external switching-off, the run-on time is also shown in the display. When this is active, it is not possible to access the KÜHLMOBIL (Mobile Chiller) via the display.

The display only then becomes re-usable when the run-on time has expired.

6.10.4 Run-on time

The run-on time is also shown on the display after switching-off. When this is active, it is not possible to access the KÜHLMOBIL (Mobile Chiller) via the display.

The display only then becomes re-usable when the run-on time has expired.

6.10.5 Flow-Meter

Is this KÜHLMOBIL equipped with a flow meter (optional), the flow quantity is shown in the display in l/ min.

The display is enabled by installing the flowmeter.

7 Maintenance

7.1 Miscellaneous

The *Maintenance* chapter comprises the areas of Care, Visual Inspection and servicing by operating personnel.

Grouping these areas into different maintenance intervals is intended to make it easier for you to plan the respectively required maintenance measures.

The instructions described in this chapter are to be understood as *minimum recommendations*.

Depending on operating conditions, expansions may be necessary in order to retain the production quality of the KÜHLMOBIL.

The indicated time intervals refer to continuous operation of the KÜHLMOBIL.



Warning!

Possible danger to persons and property assets!

As a result of incorrect inspection or servicing, direct or consequent personal injury and damage to assets are possible.

All maintenance work on the KÜHLMOBIL may only be carried out by qualified specialists, while adhering particularly to Chapter 2, *Safety*.



Note!

- Only use spare parts approved by **Van der Heijden Labortechnik GmbH**.
- **Van der Heijden Labortechnik GmbH**.

will not accept any liability if non-approved spare or replacement parts and operational equipment are used!

It must be ensured that operational materials and replaced parts are disposed of in an environmentally conscious manner.



Note!

If it is necessary to dismantle safety equipment for the purposes of maintenance work, the safety equipment must be reattached and checked immediately after completion of the work.



Note!

Follow the safety instructions and the safety instructions *for Maintenance* in Section

7.2 Maintenance and Repair

Each KÜHLMOBIL must undergo preventive maintenance measures according to the operating instructions.

The owner must ensure that the refrigeration system is checked in a satisfactory way regularly, monitored and maintained.

7.3 Operating Status

Depending on the nature and extent of the maintenance measures, the KÜHLMOBIL might have to be shut down or have its power completely disconnected.



Danger!

High voltage!

Causes death or life-threatening injuries.

The KÜHLMOBIL's mains cable absolutely must be unplugged.

7.4 Care and Visual Inspection

Task of the operator is to check the accessible areas of KÜHLMOBIL s daily for contamination and damage.

Dirt should, to the extent available, away and damage the responsible maintenance personnel or the Van der Heijden - Labortechnik GmbH are reported.

- Check the KÜHLMOBIL therefore always to starting work (see also the table visual inspection on the following page).



NOTE!

By unavoidable dirt labeling of controls or warnings may be illegible.

This can lead to incorrect operations, which can cause property damage and consequential loss.

When choosing the cleaning agent, make sure that no:

- Surfaces,
- Keyboards,
- Plastic or
- Seals

be attacked.

Unrestricted usable are all aqueous industrial cleaners.

- Clean therefore all the controls, displays and warning stickers once a week from dust and other contaminants by wiping with a damp cloth.

7.4.1 Visual Inspection

Type of Inspection	Activity/Components	✓
Check operating behaviour	Monitor KÜHLMOBIL for normal operating behaviour: – Running sounds, – heating, – development of odours. If the KÜHLMOBIL displays irregular operational behaviour, if necessary, shut down and notify maintenance personnel immediately.	
Check for residues	Check KÜHLMOBIL and surrounding area for residues of material and operating materials, and remove if necessary.	
Remove soiling	– Walkways – work surfaces – labelling, – condenser	
Check for wear	– Moving supply and drain pipes, – Externally recognisable seals, etc. – Replace components if necessary.	
Check for damage	– Externally recognisable damage to all components.	
Check for leaks	– Cooling water hoses, – Fittings, – Plug-in connections.	
Check of cooling liquid	– Clean of viewable dirt in the tank	
	–	

7.5 Servicing

7.5.1 General Instructions

All servicing work (if necessary) must be performed within the specified deadline and with the appropriate degree of care.

The components in the Mobile-Cooler housing do not require any servicing.

The only exception to this is the condenser.

- Should occur a damage, the KÜHLMOBIL ist be sent to. **Van der Heijden Labortechnik GmbH**.
- In special cases consultation with **Van der Heijden Labortechnk GmbH** be undertaken.

7.5.2 Setting the Bypass Regulating Valve

The cooling circuit of the KÜHLMOBIL is equipped with an automatically bypass regulating valve which is accessible at the rear of the housing (blue handwheel) between water flow and water return flow.

It can be used to regulate the water supply pressure by means of the flow rate, thereby influencing the flow delivery of the circulating pump.

The maximal water supply pressure can be read from the pressure gauge when the water pre-run is closed and can therefore be set to the desired value.

- The automatically bypass will be adjusted using the adjusting screw after the pump is switched off and the cap is loosed:
rotate anti-clockwise causes a reduction in pressure
rotate clockwise causes an increase in pressure.



Note!

Please take care the when adjusting the bypass regulation valve the water pressure will be enough.

Please take also care the there will not be adjusted to high pressure, which is not needed.

Adjust the required liter for the appliance which is to be cooled but do not adjust more.

The deeper the pressure on the manometer will be shown, at enough feed quantity, the better for the water pump.

7.5.3 Condenser (for Air Cooling)

Depending on the installation site and freedom from dust of the surrounding air, the condenser lamellae are to be checked for soiling at extended intervals – immediately in the event of dropping cooling performance - and cleaned, if necessary.

- To do this it is necessary to remove the front metal panel.
- It is necessary to subsequently blow the condenser from back to front using nitrogen or oil-free compressed air.
- Finally, do not forget to re-attach the metal panel.

5.4 Plate heat exchanger primary circuit

The primary circuit is at greater intervals - once with decreasing power - be checked for dirt and cleaned, if necessary. Is it not sufficient clean only the filter* to the water inlet, the internal circuit must be rinsed:

- For this purpose must be separated the primary supply hoses (water inlet / water outlet).
- Then must the primary cycle, with detergent diluted in hot water, with an external circulation pump for about 30 minutes be flushed.
- The primary supply hoses mounted (water inlet / water outlet).

Finally, don't forget the installation of the cladding sheet.



In the event of a leak in the cooling circuit can occur through the primary cooling water to an uncontrolled water outlet.

For this reason, the corresponding area on KÜHLMOBIL, in the manifold and on the device to be cooled must be monitored with moisture sensors.

In order for a leak to avoid further damage must be ensured that all water supplies and the KÜHLMOBIL immediately switched off and on again-turn, be secured.

7.5.5 Tightness test according to EN 378

"Checked for leaks" means that the device is tested by direct or indirect measuring methods primarily for leaks, in particular the parts of the device on which is most likely to leak.

Plants with a refrigerant capacity of more than 3.0kg or 1.3kg R404a R134a should be subjected to a test for leaks at least once a year.

Arise during the examination of the suspected leak, the leak with a suitable detector must be located and should be repaired and then re-examined according to national rules.

To check the tightness of several methods are used, for example, Pressurizing with inert gas (nitrogen).

Compounds have a detector or a bubble test with a detection-sensitivity in accordance with EN 1779 at max. Operating pressure, to be checked.

After Regulation 842 of the European Parliament and of the EN 378 the tightness must be checked within one month after the repair of a leak, to ensure that the repair has been effective.

7.5.6 Coolant secondary circuit

The quality of the cooling liquid in the tank from the KÜHLMOBIL shall be reviewed every 4 weeks.

Chemical or biological contaminants must be removed immediately.

This requires the container cover by unscrewing the 4 screws are removed.

The coolant circuit is then flush by the hose connections and wipe the container to clean.

The contaminated coolant must be replaced.

8 Disposal

Operation of the KÜHLMOBIL results in the accumulation of waste materials and replaced parts, which have to be disposed of in accordance with legal regulations.

8.1 Environmental Protection



Caution!

In the case of all work on the KÜHLMOBIL all legal obligations relating to waste prevention and correct recycling/disposal must be adhered too.

Especially in the case of repair and servicing work, neither the soil nor the drains must be contaminated with materials hazardous to water, such as:

- Lubricating oils and greases and
- Cleaning fluids containing solvents!

These materials must be stored, transported, collected and disposed of in suitable containers!

8.2 Coolant



Caution!

Destruction of the coolant circuit as a result of leaking coolant.

It is prohibited to open the coolant circuit without a reason and without workshop equipment suitable for the job! Opening the coolant circuit may result in it being subsequently impossible to operate the KÜHLMOBIL again.

- Work on the coolant circuit of the KÜHLMOBIL may only be performed by **Van der Heijden Labortechnik GmbH**.
- Work may only be performed by expert personnel in exceptional circumstances if there is a compelling reason for it.
- Always inform **Van der Heijden Labortechnik GmbH** before hand!
- The coolant must be correctly soaked up, and must not be allowed to escape into the environment!
- The safety specifications sheet for the coolant must be adhered to!

8.3 Final Decommissioning

If the KÜHLMOBIL is being shut down for the last time, the components and operating materials are to be disposed of in accordance with laws and regulations applicable at the time of disposal.

9 Following

9.1 Technical specifications

	Manufacturer:	
		<i>Van der Heijden Labortechnik GmbH</i>
		Tramsmeiers Berg 2
		D - 32694 Dörentrup
		Fon +49 (0) 5265 9455 20 Fax +49 (0) 5265 9455 210
		E-Mail info@van-der-heijden.de Web www.van-der-heijden.de
	General data	
☒		Product: KÜHLMOBIL Model:
☒		311-B400-MK-FF-1-7-16-38-54
☒		Power control proportional
☐		Power control w/o
☒		Temperature stability +/- 0,1 K at constant load
☐		Temperature stability +/- °C hysteresis (-K)
☒		Cooling air cooled
☒	Warranty:	2 years
☒	Case:	Art 2
☒		Wide (with wheels) 680 mm
☒		High 1530 mm
☒		Depth 750 mm
☒		Case material Sheetmetal, epoxy coating
☒		Color RAL 5003 (frame) / RAL 7035

	Type-specific data:	
☒		Coolant level, Tank ca. 124,0 liter
☒		Weight ca. 195 kg
☒		Sound pressure level 63 dB (A)
		* Measurement front in 2 m distance without reflection
	Coolant:	
☐		Medium distilled water
☐		Medium demineralized water
☒		Medium HKF 10.1
☒		Medium HKF 15.1
☒		Medium distilled water with sodium bicarbonate
☐		Medium Water-glycol mixture
☐		Medium ---
☒	Secondary connections:	As 1-Way distributor
☒		Prerun ball valve with
☒		$\frac{3}{4}$ " Hose socket
☒		Return ball valve with
☒		$\frac{3}{4}$ " Hose socket
☐	Primary connections:	
☐		Water inlet --- ---
☐		---
☐		Water outlet --- ---
☐		---
☒	Refrigerant:	R134a
☒		Capacity Approximately 3,0 kg

☒	Pressostat:	
☒		High pressure fault $\geq 18,0$ bar
☒		Low pressure fault $\leq 0,5$ bar
	Environmental conditions:	
☒		Ambient temperature, perm. +5 °C bis +32 °C
☒		Relative humidity max. 85 %
☐	Performance:	
☐	Voltage switching:	
☐		Voltage S3 = 1 ~ 110 VAC / N / PE
☐		Voltage S3 = 0 ~ 230 VAC / N / PE
☒	Fixed voltage:	
☒		Voltage ~ 400 V / 3PH / N / PE
☒		Frequency 50 Hz
☒		Amperage 14,0 A
☐		
☒		Power input 5,0 KW
☒		Cooling capacity 4.300 watts at +20 °C Water inlet temperature
☒		Operating pressure, max. 18 bar
	Circuit 1:	
☒		Pump capacity, max. 6,3 bar
☒		Pump capacity, max. 55 l/min
	Circuit 2:	
☐		Pump capacity, max. --- bar
☐		Pump capacity, max. --- l/min

	Fan:	
☒		Volume flow rate 1836 m³/h
	Thermostat parameters:	
☒		Power control: proportional
☐		Power control: thermostatically controlled tank heating
☐		Power control: w/o
☒		Voltage ~ 230 V / 50 Hz
☒		Setpoint (factory setting) +20 °C
☒		Setting range +10 °C bis +25 °C
☒		Lower deactivation 2 K under Setting range
☒		High temperature alarm +28 °C
☒		Low temperature alarm +7 °C
	Primärdaten:	
☐	Water cooled device:	
☐		Water inlet temperatue +5 °C bis +25 °C
☐		Pressure loss --- bar
☐		Max. static pressure 8,0 bar
☐		Max. differential pressure 3,0 bar
☐		Min. differential pressure --- bar
☐		Required water flow see diagram
☐	Water-water device:	
☐		Water inlet temperature --- °C
☐		ΔT --- K
☐		Required water flow --- l/h
☐		Pressure loss --- bar
☐		Max. static pressure --- bar
☐		Max. differential pressure ΔP 3,0 bar

	Options:	
☐	External switch:	
☐		Rated voltage ---
☐		potential free --- switching voltage
☒	Floating contact:	
☒		Rated voltage 250 V AC/DC
☒		Rated current 10 A AC/DC
☒	Bypass:	
☐		Gate valve manually hand adjustable
☒		Automatic opening pressure 5,5 bar
☐	Remote control:	
☐		Cable length max. ---m
	Flow display:	
☐		Circulation 1 --- l/min
☐		Circulation 2 --- l/min
☐		Circulation 3 --- l/min
☐		Circulation 4 --- l/min
☐	Operating hours counter:	
☒	Flow switch:	
☐	as pump protection:	deactivation value --- l/min
☒	as ice protection:	deactivation value 5,0 l/min
☐	In forward flow:	deactivation value --- l/min
☐	Im return flow:	deactivation value --- l/min

☒	Flowmeter:	
☐	For monitoring:	
☐	As pump protection:	deactivation value --- l/min
☐	As ice protection:	deactivation value --- l/min
☐	In forward flow:	deactivation value --- l/min
☐	In return flow:	deactivation value --- l/min
☒	To display:	
☐	In water forward flow:	deactivation value w/o
☒	In water return flow:	deactivation value w/o
☐	The total delivery rate:	deactivation value w/o
☐	For monitoring + display:	
☐	In water forward flow:	deactivation value --- l/min
☐	In water return flow:	deactivation value --- l/min
☐	The total delivery rate:	deactivation value --- l/min
☐		
☒	Time relays for overrun:	
☐	After shutdown:	Overrun time --- minutes
☐	After failure and shutdown:	Overrun time --- minutes
☒	Fan overrun	Overrun time 1 minutes
☐	Pump overrun after switching off	Overrun time --- minutes
☐	Pump overrun after fault and shutdown:	Overrun time --- minutes
☐	Overflow:	
☐	Connection:	Hose fitting ---
☐		Male thread ---
☐		With ball valve ---
☐	Automatic refill:	
☐	Connection:	Male thread ---
☐		Hose fitting ---


☒	Pressure gauge:	Anzeige	0 – 10 bar
☐	Interface RS232:	Interface	
☐	Dirt filter:		
☐		Primary	in water inlet
☐		Sekundary	in forward flow
☒	Finefilter:		
☐		Primär	in water inlet
☒		Sekundär	in water forward flow
☐	Geodetic height difference:	With solenoid valve and non-return valve	
☐	Emergency cooling:	Automatic switching to water line in case of failure	
☐	Drain:		
☐		Connection	---
☐	Drainage:	With ball valve	
☐	Internal connector:	hose	---
☐	External connector:	Hose fitting	---
☐		Male thread	---
☒	Max. fluid level		
☐		Tank to 30 liter	2,5 cm under tank rim
☒		Tank from 50liter	7,0 cm under tank rim
☒	Packing:		
☒		Material	Wooden crate
☐		Material	Palette
☐		Material	Foil
☐		Material	Box

Attention:

- The secondary circuit is isolated.
- Is there a shut down by the internal flow switch, no error will be displayed!!

9.2 Primary water consumption

(only for Chillers with Article-nr. beginning with 2-)

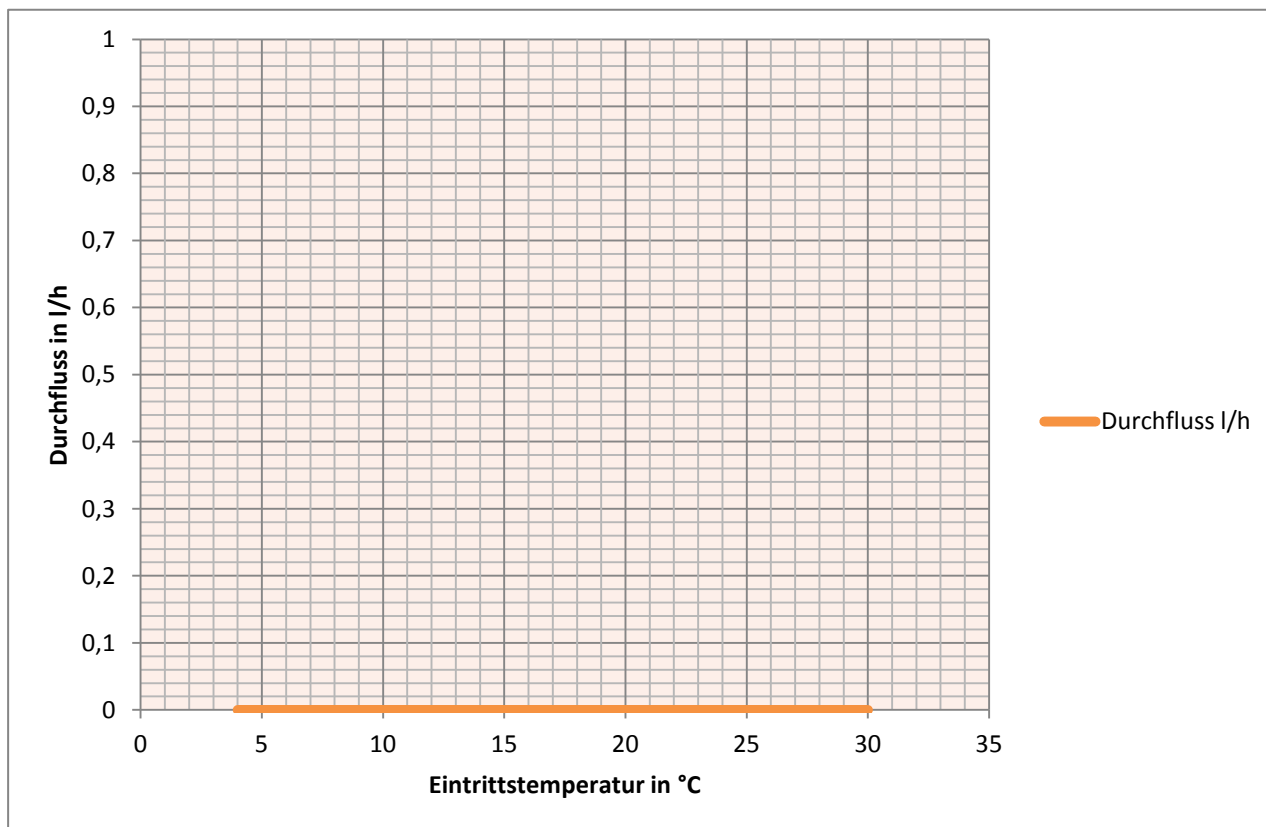

9.3 Components and spareparts

Pos.	Number	Descriptions	Article-No.	Typ	Description by wiring diagram
01	1	Fuse oil sump heater	03197	T 2,0 A	F2
02	1	Fuse terminal	16189	UK 5 - HESI	F2
03	1	Fuse	16186	T 10,0 A	F3
04	1	Fuse terminal	16190	UK 6,3 - HESI	F3
05	1	Fuse	03199	6A T	F10
06	1	Compressor	18676	MC-D8-ZB19KE-TFD	M1
07	1	Fan	---	Modell 121	M2
08	1	Pump	08525	EY-4281-MK.M5549	M3
09	1	Magnetic valve	02808	6223-A08	MV1
10	1	High pressure switch	--	see aggregat	P1
11	1	Level switch	14001	OLH-3	P2
12	1	Flow switch	13920	VH-315	P3
13	1	Flow sensor	14117	FHKU 100	P4
14	1	Low pressure switch	--	see aggregat	P5
15	1	Temperature sensor	09212	PT100	PT100

Pos.	Number	Descriptions	Article-No.	Typ	Description by wiring diagram
16	1	Main switch	15840	P1-32	Q0
17	2	Contactor	16140	DILM 12 - 10	Q1.1/Q3.1
18	1	Auxiliary switch block	16176	DILA-XHI22	Q1.1
19	2	Motor protection switch	16204	PKZM0 – 10	Q1.2/Q3.2
20	2	Standard auxiliary switch	16206	MHI-E-10-PKZ0	Q1.2/Q3.2
21	1	Oil sump heater	--	65 Watt	R1
22	1	Control unit	11501	VDH-R201	S1
23	1	Control panel	11507	VDH-B400	S2
24	1	Control unit	02820	8605	U1
25	1	Power supply	30496	UNO-PS/1AC/24DC	V1
26	1	Implement connector	14280	STASAP 2	X2
27	1	Coupling box	14300	STAK 20	X2
28	1	Filter cartridge	14656	WFMB100	---