



MCO-170AIC-PE

IncuSafe

CO₂ Incubators

165 L



Optimising cell culture outcomes and reproducibility

IncuSafe CO₂ Incubators provide precise control of CO₂ concentration and accurate, uniform, and highly responsive temperature control within the chamber. During cell culturing the inCu-saFe germicidal interior and SafeCell UV lamp work continuously to prevent contamination.

Precise & Regulated Environment


InCu-saFe and SafeCell UV both function to prevent contamination. The Direct Heat and Air Jacket System regulates the temperature whilst the Dual IR sensor controls the CO₂ level.

Time-Saving Decontamination

The high-speed decontamination system uses vaporised hydrogen peroxide and UV light. It cleans the chamber of the incubator safely in less than three hours, achieving a minimal 6 log reduction of major contaminants.

Ease of Use & Maintenance

A full colour LCD touch-screen allows full control even with gloved hands. Transfer of data is easy via a USB port. The easy-to-clean incubator interior features fully rounded corners and integrated shelf supports.

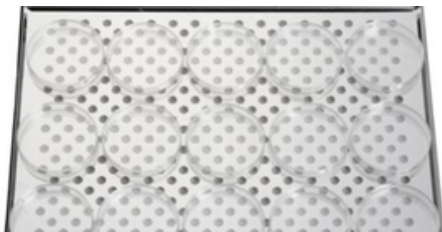

Optimum Cell Growth

Outstanding quality and performance for successful cell growth, optimal results and reproducibility. Perfect fit for the strictest and most sensitive protocols.



Efficient Workflows

Complete laboratory procedures and experiments more efficiently with less incubator downtime. Ideal for commercial applications.



Intuitive Usability

Control and visibility of the internal conditions, such as CO₂ level and temperature, is easy with the MCO-170AIC CO₂ incubator.

IncuSafe CO₂ Incubators



Direct Heat and Air Jacket System

Achieves accurate, uniform, and highly responsive temperature control within the chamber, providing exceptional uniformity and rapid recovery after door-openings.

Dual IR CO₂ Sensor


The incubator's Dual IR sensor and P.I.D control enables ultra-fast CO₂ recovery without overshoot, even following multiple door-openings.

Active Background Decontamination

The exclusive InCu-saFe copper-enriched stainless steel alloy interior offers the germicidal properties of copper and the durability of stainless steel. The optional, isolated, SafeCell UV lamp decontaminates circulating air and water in the humidifying pan, without harming cultured cells.

Condensation Management

With a unique antibacterial coating, the 'dew stick'—controlled by Peltier technology—condenses water on its surface, which then drips into the humidifying pan, preventing unwanted condensation in the chamber and possible contamination.

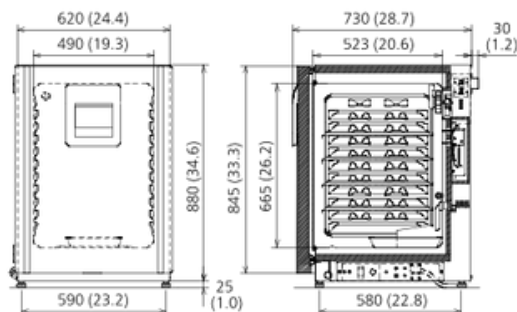

The MCO-170AIC series is certified as a Class IIa Medical Device (93/42/EEC and 2007/47/EC) for medical purposes of culturing cells, tissues, organs and embryos.



PHC Europe B.V.
Nijverheidsweg 120 | 4879 AZ Etten-Leur | Netherlands
T: +31 (0) 76 543 3839 | F: +31 (0) 76 541 3732
www.phc.com/eu/biomedica

Model Number		MCO-170AIC-PE	MCO-170AICUV-PE	MCO-170AICUVH-PE
External Dimensions (W x D x H) ¹⁾	mm	620 x 730 x 905		
Internal Dimensions (W x D x H)	mm	490 x 523 x 665		
Volume	liters	165		
Net Weight	kg	80		
Performance				
Temperature Control Range & Fluctuation	°C	AT +5 ~ +50, ±0.1		
Temperature Uniformity ²⁾	°C	±0.25		
CO ₂ Control Range & Fluctuation	%	0 ~ 20, ±0.15		
Humidity Level & Fluctuation	%RH	95, ±5		
Sterilisation Method		H ₂ O ₂ Decontamination		
Control				
Temperature Sensor		Thermistor		
CO ₂ Sensor		Dual IR		
Display		LCD Touch Screen		
Construction				
Exterior Material		Painted Steel (rear cover not painted)		
Interior Material		Stainless Steel Copper-Enriched Alloy		
Insulation Material		Extruded polystyrene		
Heating Method		Direct Heat & Air Jacket System		
Outer Door		1		
Outer Door Lock	qty	1		
Field Reversible Door		Optional	Optional	Standard
Inner Door		Included		
Shelves	qty	1 gastight - made of tempered glass		
Shelf Dimensions (W x D x H)	qty	4 x Stainless Steel Copper-enriched Alloy		
Max. Load per Shelf	mm	470 x 450 x 12		
Max. Shelf Capacity	kg	7		
Access Port	qty	10		
Access Port Position	qty	1		
Access Port Diameter		Rear Upper Left		
Alarms	Ø mm	30		
(R = Remote Alarm, V = Visual Alarm, B = Buzzer Alarm)				
Power Failure				
Out of Temperature Setting		R		
High Temperature		V-B-R		
Out of CO ₂ Setting		V-B-R		
Door open		V-B-R		
Electrical and Noise Level		V-B		
Power Supply				
Frequency	V	230		
Noise Level ³⁾	Hz	50		
Options	dB	29		
SafeCell UV® System				
		MCO-170UVS-PE ⁴⁾	Standard	
H ₂ O ₂ Decontamination Board		MCO-170HB-PE ⁴⁾		Standard
Electric Door Lock with Password		MCO-170EL-PW ⁴⁾		Standard
H ₂ O ₂ Vapor Generator		MCO-HP-PW ⁴⁾		
H ₂ O ₂ Reagent, pack of 6 bottles		MCO-H202-PE		
Multiple Inner Doors		MCO-170ID-PW		
CO ₂ Gas Pressure Regulator		MCO-100L-PW		
Automatic CO ₂ Cylinder Changeover System		MCO-21GC-PW		
Semi-automatic one point Gas Calibration Kit		MCO-SG-PW		
InCu-saFe® Shelf		MCO-170ST-PW		
InCu-saFe® Half Tray System		MCO-25ST-PW		
Double Stacking Bracket*		MCO-170PS-PW		
Stacking Plate*		MCO-170SB-PW		
Roller Base		MCO-170RB-PW		
Optional communication systems ⁵⁾				
Ethernet interface (LAN)		MTR-L03-PW		
Digital interface (RS232C/RS485)		MTR-480-PW		

¹⁾ Exterior dimensions of main cabinet only, excluding handle and other external projections.

²⁾ ±0.25°C; ambient temp 23°C - 25°C, setting 37°C, CO₂ 5%, no load

³⁾ Nominal value

⁴⁾ MCO-170AIC series requires MCO-170HB-PE, MCO-170EL-PW, MCO-HP-PW and SafeCell UV option for H₂O₂ decontamination

⁵⁾ MCO-170AIC series can only be fitted with one communications interface.

* If stacking two incubators, make sure the double-stacking dedicated securing hardware and spacer are used.