

## Data sheet

## Laboratory Sterilizer SRW 1000 Smart



The photo above is for reference only, may show additional options not included in standard equipment. The real appearance, particularly color and structure of the material may differ from the ones presented in the photo.

### Advantages of the SMART controller:

- 4,3", clear, full colour touch screen
- LAN, USB ports for data transfer
- multi-segment time and temperature programs
- visual and sound alarm
- internal memory for programs and data storage
- event registry
- user manual for direct download
- Quick change of program parameters
- Alarm Bar
- operating with gloves on



**Smart** - preview screen

## TECHNICAL DATA

|                      |                               |
|----------------------|-------------------------------|
| air convection       | forced                        |
| chamber capacity [l] | 1005                          |
| working capacity [l] | 1005                          |
| controller           | microprocessor PID            |
| display              | 4,3" full colour touch screen |

## TEMPERATURE


|  |   |
|--|---|
| temperature range [°C]                   | 5°C above ambient temperature ... +250°C    |
| temperature resolution every ... [°C]    | 0,1   |
| temperature fluctuation at 105°C [+/°C]* | 0,6   |
| temperature variation at 105°C [+/°C]*   | 3,0   |
| over temperature protection              | class 2.0 to DIN 12880 / class 3.1 (option) |

## CHAMBER

|                 |   |
|-----------------|---|
| door            | solid / door with viewing window (option) |
| <b>interior</b> |   |
| Smart           | acid-proof stainless steel to DIN 1.4301  |
| IG Smart        | acid-proof stainless steel to DIN 1.4301  |
| <b>housing</b>  |   |
| Smart           | powder coated sheet                       |
| IG Smart        | stainless steel linen finish              |

## overall dims [mm] /1/

|                           |      |
|---------------------------|------|
| A width                   | 1260 |
| B height                  | 2000 |
| depth C                   | 880  |
| <b>internal dims [mm]</b> |      |
| width D                   | 1040 |
| height E                  | 1610 |
| depth F                   | 600  |

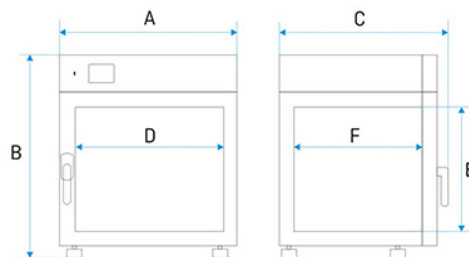

|  |        |
|--|--------|
| shelves (standard   max)                 | 6   22 |
| - reinforced shelf version (PW) [kg] /3/ | 100    |
| max unit workload [kg]                   | 300    |
| weight [kg]                              | 370    |

## ELECTRICAL PARAMETERS

|                   |                   |
|-------------------|-------------------|
| power supply      | 400V 50-60Hz      |
| nominal power [W] | 5500              |
| warranty          | 24 months         |
| manufacturer      | POL-EKO-APARATURA |

all the above technical data refer to standard units (without optional accessories)

\* - fluctuation measured in centre of the chamber; in space, variation (K) calculated for chamber as:

$K = \frac{T_{\text{average max.}} - T_{\text{average min.}}}{2}$

\*\* - other power supplies on request


1 - depth doesn't include 50 mm of power cable

2 - on uniformly loaded surface

3 - reinforced shelf

4 - reinforced version

## OPTIONS AND ACCESSORIES

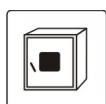


Order number: \*/A

Door with viewing window




Order number: \*/P INOX

Stainless steel wire shelf INOX

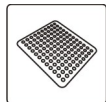

Order number: \*/PP

Perforated shelf

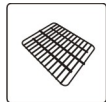

Order number: \*/PW

Reinforced shelf



Order number: KUW GN\*/\*

Stainless steel cuvettes




Order number: OWW/OWW LED

Interior lighting



Order number: HEPA

HEPA Clean Air Filter

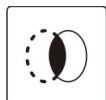


Order number: OCZ/N

Non-standard access port for external sensor



Order number: LabDesk

LabDesk software

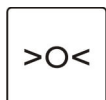

Order number: BRT\*/L or IQ/OQ/PQ

Calibration and IQ, OQ, PQ qualification



Order number: \*/3.1

Over temperature protection 3.1 class according to DIN 12880



tel. (+48) 32 453 91 70

e-mail: [info@pol-eko.com.pl](mailto:info@pol-eko.com.pl)  
internet: [www.pol-eko.com.pl](http://www.pol-eko.com.pl)

## POL-EKO-APARATURA SP. J.

A. POŁOK – KOWALSKA, S. KOWALSKI  
44-300 WODZIŚLAW ŚLĄSKI, UL. KOKOSZYCKA 172 C



Order number: BPP 12

Battery backup for display



Order number: PORT ALARM

Dry alarm contact