

G 7804

Laboratory glassware washing machine



G 7804

for laboratory glassware



The G 7804 Lavabor has been designed for cleaning laboratory glassware in schools, industrial, environmental and research laboratories.

Capacity

- 39 narrow necked flasks (E 329) or
- 116 pipettes (E 406) or
- 1600 test tubes (E 103/E 139)

Construction

- H 850 (820), W 600, D 600 mm
- Front loading machine with drop down door
- Suitable for use freestanding or can be fitted under a worktop
- Stainless steel casing and lid
- Space frame construction with side wall insulation
 - low heat and noise emissions
 - easy to recycle
- Wash chamber and spray arms in high grade stainless steel

Features

- Wash chamber with 2 levels
 - high capacity
- Two spray arms (3rd spray arm is in the upper basket)
 - accurate wash
- Direct coupling of upper basket and injector mobile units to water feed
- Profi Monobloc Water softener
 - continual reactivation during the cleaning process
- Water intake controlled by fly wheel counter
 - precise amount of water taken in ensures the correct ratio of water to cleaning and disinfecting agents
 - automatic extension in low pressure installations
 - booster pump not necessary
 - consumption linked reactivation of water softener
- Electro thermal door lock
 - safety for personnel
 - high process safety
- Access to check and control temperatures and temperature holding times

Dispenser system

- One dispenser for powder cleaning agents and two for liquid agents in the door

Options and accessories

- External pump for non-pressurised demineralised water feed
- Dispenser module (DOS G60 or DOS K60) for liquid cleaning agents
- Dispenser (DOS G10) for neutralising liquids
- Plinth H 300 mm



G 7804 controls

Controls

Multitronic Novo Vario

Programmes

- 3 standard cleaning programmes

Programme Control

- Single switch control knob

Display

- Programme sequence display
- Temperature and programme duration
- Programme end, optical and acoustic signal
- Service and fault check lights

Interface

- Optical interface for service technician



Programme chart

Programme	When to use	Pre-wash
A SHORT	For laboratory utensils and glassware with light water soluble contamination. This programme does not include a pre-wash or neutralisation. The final rinse is carried out with de-ionised (AD) water.	
B MEDIUM	For normally soiled laboratory utensils and glassware where the contamination is soluble in water or at higher temperatures where alkaline. This programme does not include a pre-wash. The final rinse is carried out with de-ionised (AD) water.	
C LONG	For heavily soiled laboratory utensils and glassware where the contamination is soluble in water or at higher temperatures where alkaline. This programme includes a pre-wash and neutralisation. The final rinse is carried out with de-ionised (AD) water.	X
DRAIN	For draining water out of the machine, e.g. when a programme has been interrupted. Turn the programme selector to the STOP position (⏏) first!	
PRE-RINSE	For rinsing heavily soiled items (e.g. for removing soiling or disinfectant residue or to prevent soiling drying on to glassware when a complete programme does not yet need to be run).	

X = Sections included in a programme

¹⁾ Temperatures/temperature holding times in the main wash and final rinse stages can be used to suit user requirements

²⁾ Rinse 2 has to be activated by the user if required. The machine can be programmed by a service technician to use AD water for this rinse if required.

DOS 1 = For dispensing cleaning agent

DOS 2 = For dispensing rinsing agent

DOS 3 = For dispensing neutralising agent

AD = Aqua destillata (programmable for Rinse stage ²⁾)



	Main wash ¹⁾	Rinse 1	Rinse 2 ²⁾	Final rinse ¹⁾	Drying (optional)
	X DOS 1 70°C / 5'		(X) DOS 3 (AD)	X 70°C / 3' AD (DOS 2)	(X)
	X DOS 1 70°C / 5'	X DOS 3	(X) (AD)	X 70°C / 3' AD (DOS 2)	(X)
	X DOS 1 70°C / 5'	X DOS 3	(X) (AD)	X 70°C / 3' AD (DOS 2)	(X)
		X			

Baskets and inserts

Superb cleaning results using special baskets and inserts

The Miele fresh water circulation system, coupled with the double wide-jet spray arms washes on two levels.

The upper and lower stainless steel baskets, O 188 and U 874 can be fitted with a wide range of inserts to suit individual wash requirements.




O 188/2 upper basket (empty)

- Open at front
- For various inserts
- H 215, W 531, D 475 mm
- Clearance: 165 mm
- Height adjustable +/- 20 mm
- Built-in spray arm


O 190/2 upper basket (empty)

- As O 188/2
- Clearance: 215 mm

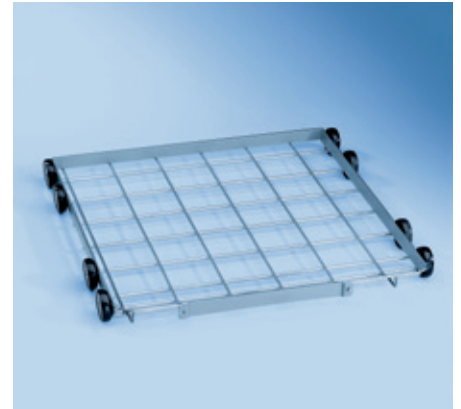

U 874/1 lower basket (empty)

- Open at front
- For various inserts
- H 50, W 534, D 515 mm
- Clearance: 270 mm with O 188/2
- Clearance: 220 mm with O 190/2
- Height adjustable +/- 20 mm

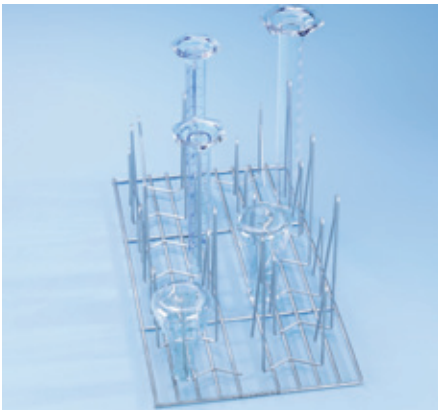

E 106 Insert 1/2

- With 28 large and small spring clips
- For lower basket

E 106/1 Insert 1/2

- With 28 x 105 mm small spring clips
- For upper or lower basket

E 106/2 Insert 1/2

- With 15 x 175 mm large spring clips
- For lower basket



E 109 Insert 1/2

- For 21 beakers up to 250 ml, round and conical flasks etc.
- For lower basket

E 110 Insert 1/2

- For 10 beakers 250–600 ml.
- For lower basket

E 111 Insert 1/2

- For 8 beakers 600–1000 ml.
- For lower basket

E 144 Insert 1/2

- For 18 beakers up to 250 ml.
- For upper or lower basket



E 118 Insert 1/1

- For 38 Petri dishes
- For upper or lower basket

E 136 Insert 1/1

- For 56 Petri dishes
- For lower basket

E 137 Insert 1/1

- For 56 Petri dishes
- For use with E 136

Baskets and inserts



Segment inserts

103/1 Insert 1/4

- For approx. 200 x 75 mm long test tubes
- For upper or lower basket

E 104/1 Insert 1/4


- For approx. 200 x 105 mm long test tubes
- For upper or lower basket

E 105/1 Insert 1/4

- For approx. 200 x 165 mm long test tubes
- For lower basket

E 139/1 Insert 1/4

- For approx. 200 x 200 mm long test tubes
- For lower basket

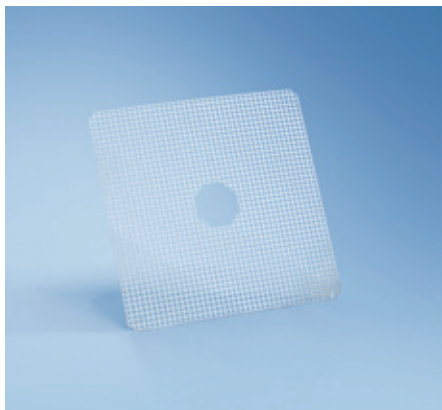

A 11 Insert 1/1

- Stainless steel mesh underlay
450 x 450 mm
- For upper or lower basket

A 12 Insert 1/2

- Stainless steel mesh underlay
450 x 225 mm
- For upper or lower basket



A 13 Lid

- For inserts E 103/1, E 104/1, E 105/1 and E 139/1
- Stainless steel



E 149 Insert 1/4

- For approx. 80 x 105 mm long test tubes
- For upper or lower basket



AK 12 Insert 1/2


- Mesh basket for funnels, beakers, wide-necked vessels, etc.
- For upper or lower basket



A 14 Lid


- For insert AK 12
- 210 x 210 mm
- Stainless steel

Baskets and inserts

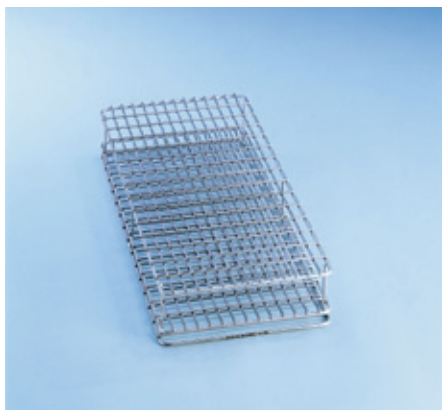

E 134 Insert 1/2

- For 210 slides
- For upper or lower basket

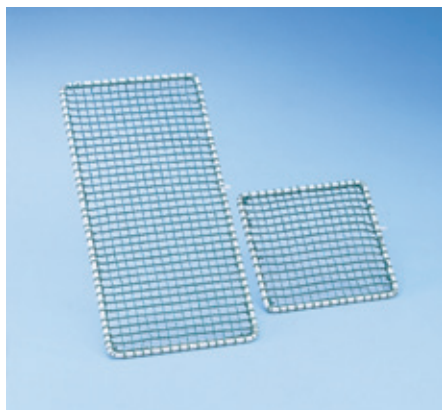

A2 Cover net 1/2

- 216 x 456 mm plastic coated metal frame with plastic netting for half inserts

A3 Cover net 1/4

- 206 x 206 mm plastic coated metal frame with plastic netting for quarter inserts



E 403 Insert 1/2

- For 105 Petri dishes, 50–60 mm Ø
- For upper or lower basket

E 402 Insert 1/2

- For 44 Petri dishes, 80–125 mm Ø
- For upper or lower basket



E 329 Mobile injector unit 1/1

- For narrow-necked glassware with:
4/5 jets 2.5 mm Ø x 90/110 mm
5/5/5 jets 4 mm Ø x 140/160/180 mm
5/5/5 jets 6 mm Ø x 200/220/240 mm
1 jet for flushing out powder dispenser
- For lower basket



E 340 Mobile injector unit

Multi-purpose

- Half for narrow-necked glassware with:
3/3/3 jets 4 mm Ø x 140/160/180 mm
3/3/4 jets 6 mm Ø x 200/220/240 mm
1 jet for flushing out powder dispenser
- Half for other inserts
- For lower basket



E 350 Mobile injector unit 1/1

- For narrow-necked glassware with:
15 jets 4 mm Ø x 160 mm
15 clips
18 jets 6 mm Ø x 220 mm
18 clips
1 jet for flushing out powder dispenser
- For lower basket

Baskets and inserts


Accessories



E 355 Mobile injector unit

Multi-purpose

- Half for narrow-necked glassware with:
7 jets 4 mm Ø x 160 mm
7 clips
9 jets 6 mm Ø x 220 mm
9 clips
1 jet for flushing out powder dispenser
- Half for other inserts
- For lower basket

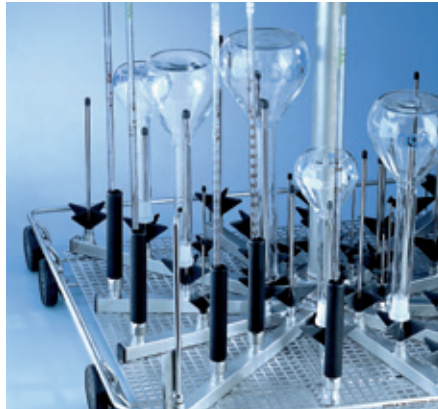


E 336 MIBO pipette holders

- Plastic with screw fitting for holding individual pipettes in a mobile injector unit



E 406 Mobile injector unit

- For 116 pipettes up to 450 mm
- Complete with support frame

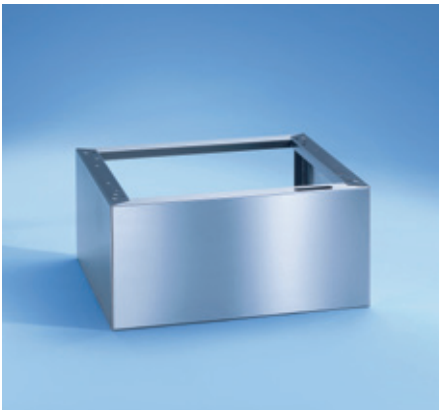

UE 30-60/60-78 plinth

- H 300 x W 600 x D 600 mm
- Stainless steel
- Can be bolted to the machine




G 60 Liquid detergent dispensing pump

- Shown above with container (not supplied)

Technical data G 7804 Lavabor

Controls

Multitronic Novo Plus	•
Impulse heating for drying, T = 10 min.	•
Water connections	
1 x cold / hot, 0.5–10 bar flow pressure (50–1000 kPa)	•
1 x AD (VE) water, 0.5–10 bar flow pressure (50–1000 kPa)	•
2 inlet hoses 1/2" with 3/4" threaded union, length approx. 1.7 m	•
1 x drain hose, length approx. 1.4 m, internal dia. 22 mm – delivery head: 100 cm	•
Circulation pump Qmax (l/min)	400
Electrical connection	
1 N AC 230 V 50 Hz	1•
Heating/Wash cabinet	2.5
Circulation pump [kW]	0.7
Total rated load [kW]	3.2
Fuse rating (AUS version) [A]	1 x 15
Mains connection cable, length approx. 2.0 m, Ø 3 x 1.5 mm ²	•
Dispenser system	
1 powder dispenser in the door	•
2 x liquid dispensers in the door (switched on together)	•
Connections for	
DOS G 10 dispenser for liquid rinsing agents	•
DOS G 60 / DOS K60 dispenser for liquid cleaning agent	•
Water softener	
For cold or hot water to max 70°C, Monobloc	•
Dimensions, weight	
External dimensions: H 850 (820), W 598, D 600 mm	•
Weight, unloaded [kg]	70

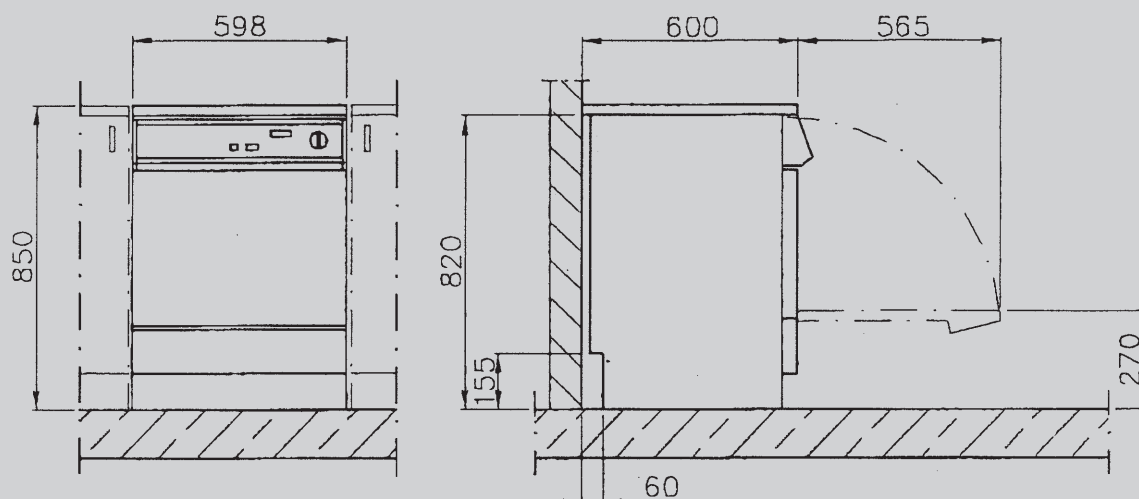

Other Miele Professional products for different target groups



PW 6065 Vario washer extractor

- 6.5 kg load capacity
- Sluice and thermal disinfection programme
- Dryer can be stacked on top
- Maintenance free motor
- 1400 rpm spin speed



For bigger applications:

Where the throughput demands a greater load capacity than the Little Giants, Miele has the answer. The capacity range extends to 32 kg and even includes barrier washing machines for specialist applications.



Washer disinfectors for dental surgeries G 7881

- Cleans and disinfects
- External dimensions:
H 850 (820), W 600, D 600 mm
- Electrical connection:
3 N AC 400 V 50 Hz
- Electronic controls with programme buttons for 5 programmes



For hospitals:

PG8536

Washer disinfector

- Control system: Profitronic+, freely programmable
- Integrated hot-air drying unit and storage drawer for liquid agent supply containers
- Capacity: 3 AN sets or 2 MIS sets or 7 DIN mesh trays
- H 1175, W 900, D 700 mm



G 7826 Washer disinfector

- Control system: Freely programmable "Profitronic" controls
- Two door barrier machine
- Wash cabinet with drop down door
- Batch capacity:
10 DIN mesh trays
or 4 AN sets
or 2 MIS sets
- External dimensions (plinth, washer disinfector, drying unit):
H 1975, W 900, D 750 mm

PG 8528 Washer disinfector

- Control system: Freely programmable "Profitronic" controls
- Two door barrier machine
- Wash cabinet with lift-up door
- Batch capacity:
15 DIN mesh trays
or 7 AN sets
or 3 MIS sets
- External dimensions (plinth, washer disinfector, drying unit, and if applicable, steam condenser):
H 2420, W 1150, D 870 mm

Miele Australia Pty Ltd
1 Gilbert Park Drive
Knoxfield VIC 3180

Tel: 1300 464 353
Fax: (03) 9764 7149
www.miele-professional.com.au
info@miele-professional.com.au

Miele & Cie. KG
Carl-Miele-Straße 29, 33332 Gütersloh
Postfach, 33325 Gütersloh

www.miele-professional.com
info@miele.de

Miele worldwide:

AUSTRALIA
Miele Australia Pty. Ltd.
1 Gilbert Park Drive
Knoxfield, VIC, 3180
Tel. 1300 464 353
Fax (03) 9764 7149

AUSTRIA
Miele Gesellschaft mbH.
Mielestraße 1
5071 Wals/Salzburg
Tel. +43 50 800-81 481
Telefax +43 50 800-81 429

BELGIUM
N.V. Miele S.A.
z.5 Mollem 480
Hof te Bollebeeklaan 9
1730 Mollem
Tel. +32 2451 14 11
Telefax +32 2451 14 14

CANADA
Miele Limited
161 Four Valley Drive
Vaughan, ON L4K 4V8
Tel. +1 905 532 2270
Telefax +1 904 532-2290

CHILE
Miele Electrodomesticos Ltda.
Nueva Costanera 4055
Vitacura
Santiago de Chile
Tel. +56 2 95700-00
Telefax +56 2 95700-79

CHINA
Miele (Shanghai) Trading Ltd
Beijing Branch
3F South Tower, Kerry Center
1 Guang Hua Road
Chaoyang District
100020 Beijing
Tel. +86 10 5783 2688
Telefax +86 10 5783 2600

Miele (Shanghai) Trading Ltd.
1-3 Floor, No. 82 Shi Men Yi Road,
Jing An District
200040 Shanghai PRC
Tel. +86 21 6135 3559
Telefax +86 21 6135 3993

CROATIA
Miele trgovina i servis d.o.o.
Buzinski prilaz 32
10000 Zagreb
Tel. +385/1/66 89 000
Telefax +385/1/66 89 090

CZECH REPUBLIC
Miele spol. s.r.o.
Holandská 4
63900 Brno
Tel. +420 5 4 35 53-1 11
Telefax +420 5 4 35 53-1 19

DENMARK
Miele A/S
2600 Glostrup, Erhvervsvej 2
Tel. +45 43 27 11 00
Tel. Salg/Professional
+45 43 27 15 00
Tel. Jylland +45 97 12 70 66
Telefax hovednr. +45 43 27 11 09
Telefax Salg/Professional
+45 43 27 15 09

ESTONIA
Miele SIA Eesti filiaal
A.H. Tammsaare tee 47
11316 Tallinn
Tel. +371 2922 3571

FINLAND
Miele Oy
Porttikaari 6
01200 Vantaa
Tel. +3589 87 59 70
Telefax +3589 87 59 72 99

FRANCE
Miele S.A.S.
9 avenue Albert Einstein
93151-Le Blanc-Mesnil (Paris)
Tel. +33(0)1/49 39 34 44
Telefax +33(0)1/49 39 44 38

GERMANY
Miele & Cie. KG
Vertriebsgesellschaft Deutschland
Carl-Miele-Str. 29
33332 Gütersloh
Tel. 0800 22 44 644*
Telefax 0800 33 55 533*
*Only in Germany

GREAT BRITAIN
Miele Co. Ltd.
Fairacres
Marcham Road, Abingdon
Oxon OX14 1TW
Tel. 0845 330 3618
Telefax 01235 405636

GREECE
Miele Hellas E.P.E.
Mesogion 257
154 51 N. Psychiko-Athens
Tel. national 01/10 6 79 44 44
international 00 30/2/10 6 79 44 44
Telefax 01/10 6 79 42 09
international 00 30/2/10 6 79 42 09

HUNGARY
Miele Kft.
Alsó Törökvezs út 2.
1022 Budapest
Tel. +36 1 8806-400
Telefax +36 1 8006-402

HONG KONG
Miele (Hong Kong) Ltd.
41/F - 4101, Manhattan Place
23 Wang Tai Road,
Kowloon Bay
Tel. +852/2610 1331
Telefax +852/3579 1404

INDIA
Miele India Pvt. Ltd.
1st Floor,
Copia Corporate Suites,
Plot No. 9, Jasola,
New Delhi - 110025
Tel. +91 (0)11 469000-00
Telefax +91 (0)11 469000-01

IRELAND
Miele Ireland Ltd.
2024 Bianconi Avenue
Citywest Business Campus
Dublin 24
Tel. +3531 461 07 10
Telefax +3531 461 07-97

ITALY
Miele Italia S.r.l.
Strada di Circonvallazione, 27
39057 S. Michele-Appiano (BZ)
Tel. +39 04 71/66 61 11
Telefax +39 04 71/66 63 50

JAPAN
Miele Japan Corp.
Meguro Yamate Place 9F
2-10-11 Meguro
Meguro-ku, Tokyo
153-0063 Japan
Tel. +81(3)5740-0034
Telefax +81(3)5740-0035

KAZAKHSTAN
Miele LLP
Building 1B (office 602)
Al-Farabi Pr. 13
050059 Almaty
Tel. +7 727 3111141
Telefax +7 727 3111042

KOREA
Miele Korea Limited
Miele House 8 fl.,
607-10 Yeoksam-dong
Gangnam-gu
Seoul 135-080
Tel. +82 2 3451 9333
Telefax +82 2 3451 9399

LATVIA
Miele SIA
Zaubes Str. 9A-38
1013 Riga
Tel. +371 66016600
Telefax +371 67889879

LITHUANIA
SIA "Miele" Lietuvos filialas
Laisvės pr. 125
06118 Vilnius
Tel. +371 29223521

LUXEMBOURG
Miele s.à r.l.
20, Rue Christophe Plantin
2339 Luxembourg
Tel. +352/4 97 11-25
Telefax +352/4 97 11-39

MALAYSIA
Miele Sdn Bhd
Suite 12-2, Level 12
Menara Kencana Petroleum
Solaris Dutamas
No. 1, Jalan Dutamas 1
Kuala Lumpur 50480
Tel. +603 6209 0288
Telefax +603 6205 3768

MEXICO
Miele S.A. de C.V.
Av. Santa Fé 170
German Centre 0-4-2
Col. Lomas de Santa Fé
C.P. 01210 México, D.F.
Tel. +52/55 85 03 98 70
Telefax +52/55 85 03 98 74

NETHERLANDS
Miele Nederland B.V.
De Limiet 2
4131 NR Vianen
Tel. +31 3 47/37 88 83
Telefax +31 3 47/37 84 29

NEW ZEALAND
Miele New Zealand Ltd.
Unit L, 10-20 Sylvia Park Road
Mt. Wellington, 1060 Auckland
Tel. +64 9573 1269
Telefax +64 9573 1268

NORWAY
Miele AS
Nesbruveien 71
1394 NESBRU
Tel. +47/67 17 31 00
Telefax +47/67 17 31 10

POLAND
Miele Sp. z. o. o.
ul. Gotarda 9
02-683 Warszawa
Tel. +48 22/5 48 40 00
Telefax +48 22/5 48 40 10

PORTUGAL
Miele Portuguesa, Lda.
Av. do Forte, 5
2790-073 Carnaxide
Tel. +351/21/42 48-100
Telefax +351/21/42 48-109

ROMANIA
Miele Appliances SRL
Intrarea Narciselor 8
Otopeni
Jud. Ilfov
Tel. +40 21 3520777
Telefax +40 21 3520776

RUSSIA
000 Miele CIS
Leningradsky Prospekt, 31A, bld. 1
125284 Moscow
Tel. +7 495 745-8992
Telefax +7 495 745-8680

SERBIA
Miele d.o.o. Beograd
Balkanska 2
Beograd - Stari Grad 11000
Tel. +381 11 2055060
Telefax +381 11 2055065

SINGAPORE
Miele Pte. Ltd.
163 Penang Road
#04-02/03 Winsland House II
Singapore 238463
Tel. +65/67 35 11 91
Telefax +65/67 35 11 61

SLOVAKIA
Miele s. r. o.
Plynárenská 1
82109 Bratislava
Tel. +421/2/58 10 31 11
Telefax +421/2/58 10 31 19

SLOVENIA
Miele d.o.o.
Brnčičeva ulica 41 g
1231 Ljubljana - Črnuče
Tel. +386 1563 44-80
Telefax +386 1563 44-90

SOUTH AFRICA
Miele (Pty) Ltd.
63 Peter Place
Bryanston 2194
Tel. +27(0)11/875-9000
Telefax +27(0)11/875-9035

SPAIN
Miele S.A.U.
Avda. de Bruselas, 31
28108 Alcobendas (Madrid)
Tel. +34/91/6 23 20 00
Telefax +34/91/6 62 02 66

SWEDEN
Miele AB
Industrivägen 20
Box 1397
171 27 Solna
Tel. +46 85 62 29-000
Telefax +46 85 62 29-209

SWITZERLAND
Miele AG
Limmatstr. 4
8957 Spreitenbach
Tel. +4156 4 17 20 00
Telefax +4156 4 17-24 69

TURKEY
Miele Elektrikli Aletler
Diş. Tic. ve Paz. Ltd. Şti.
Güvercin Sokak No. 10
34330 Levent-Istanbul
Tel. +90/2 12/3 86-08 00
Telefax +90/2 12/3 25 84 49

UKRAINE
Miele LLC
Zhylyanskaya Str. 48, 50 A
01033 Kiev
Tel. +38 044 4960300
Telefax +38 044 4942285

UNITED ARAB EMIRATES
Miele Appliances Ltd.
Gold & Diamond Park
Sheikh Zayed Road
Building No. 6,
Offices No. 219
Dubai
Tel. +971 4 3418444
Telefax +971 4 3418852

USA
Miele, Inc.
9 Independence Way
Princeton, NJ 08540
Tel. 001-609/419 43 74
Telefax 001-609/419 18 53

OTHER COUNTRIES:
Sales International
Miele & Cie. KG
Carl-Miele-Straße 29
33332 Gütersloh
Tel. +49(0)52 41-89 15 07
Telefax +49(0)52 41-89 15 00