



Speed, Performance, Value - The Right Balance


The Adam Anti-Vibration Table helps provide a stable surface and reduces vibrations, shock and outside influences effecting devices placed on the marble working surface. Designed for the laboratory or metrology room the anti-vibration table is ergonomically styled yet built for durability.

- Large working surface of 400mm x 450mm
- Easy to set up light weight mild steel tubular frame
- Adjustable levelling feet
- Scratch resistant marble slab
- Polished finish

The table comprises two separate parts:

- An internal working surface consisting of a large marble slab mounted on adjustable rubber shock absorbing mounts.
- An external service table with steel frame and adjustable levelling feet.

Overall Dimensions (w x d x h)	800 x 600 x 782 mm
Work Surface Dimensions	400 x 450 mm
Weight of Work Surface	48 kg
Total Weight	82 kg

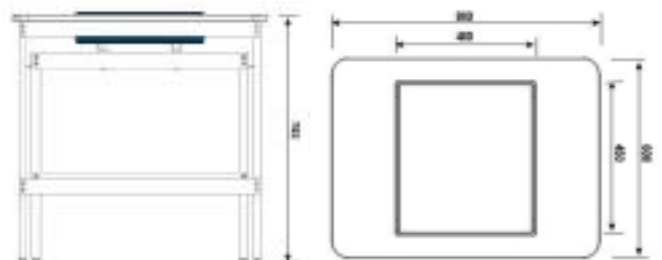

www.adamequipment.com

Adam Equipment - offices worldwide