



POL-EKO®

(+48) 32 453 91 70
info@pol-eko.com.pl
www.pol-eko.com.pl

POL-EKO A. Polok-Kowalska sp.k.
44-300 Wodzisław Śląski
ul. Kokoszycka 172 C

Data sheet

print date: 2026-07-20

Laboratory Freezer ZLN-T 125 Smart PRO



The photo above is for reference only, may show additional options not included in standard equipment. The real appearance, particularly color and structure of the material may differ from the ones presented in the photo.

Advantages of the SMART PRO controller:

- large (7"), clear, full colour touch screen
- LAN, USB ports and WiFi for communication and data transfer
- multi-segment time and temperature programs
- overview of data in tabular and graphic form
- visual and sound alarm
- Admin function for management
- password protected log-in
- internal memory for programs and data storage
- event registry with user notifications
- LabDesk software and user manual for direct download



Smart PRO - preview screen



POL-EKO®

(+48) 32 453 91 70
info@pol-eko.com.pl
www.pol-eko.com.pl

POL-EKO A.Polok-Kowalska sp.k.
44-300 Wodzisław Śląski
ul. Kokoszycka 172 C

TECHNICAL DATA

air convection	natural
chamber capacity [l]	138
working capacity [l]	127
controller	microprocessor PID
display	7" full colour touch screen
insulation type	polyurethane foam

TEMPERATURE

temperature range [°C]	-40...0
temperature resolution every ... [°C]	0,1
temperature fluctuation at -20°C [±/°C]*	0,5
temperature variation at -20°C [±/°C]*	2
temperature protection	class 3.2 to DIN 12880 (option)

CHAMBER

door type	solid
-----------	-------

interior

P Smart PRO	corrosion-resistant stainless steel to DIN 1.4301
PS Smart PRO	corrosion-resistant stainless steel to DIN 1.4301

housing


P Smart PRO	powder coated sheet
PS Smart PRO	stainless steel polished

overall dims [mm] / 1/

width A	720
height B	1220
depth C	810

internal dims [mm]

width D	390
width D'	420
height E	600
depth F	530
depth F'	550

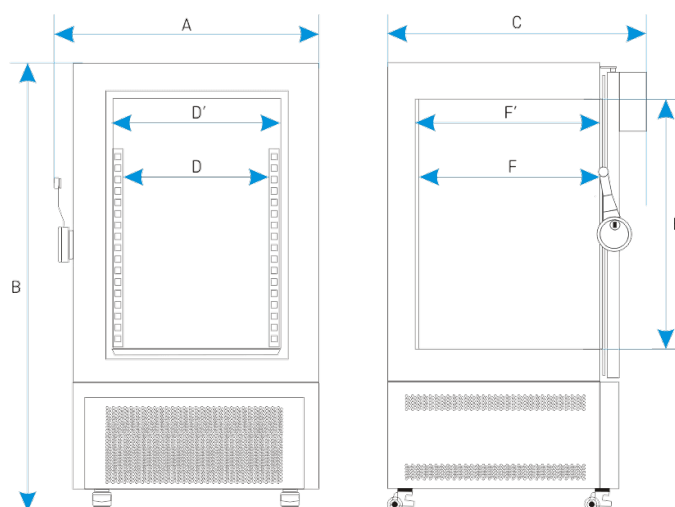

shelves (standard max)	2 3
max shelf workload [kg] /2/	10
- reinforced shelf version (PW) [kg] /3/	50
max unit workload [kg]	50
- reinforced unit version (W) [kg] /4/	100
weight [kg]	105



POL-EKO®

(+48) 32 453 91 70
info@pol-eko.com.pl
www.pol-eko.com.pl

POL-EKO A. Polok-Kowalska sp.k.
44-300 Wodzisław Śląski
ul. Kokoszycka 172 C

ELECTRICAL PARAMETERS

voltage	230V 50/60Hz
nominal power [W]	450
refrigerant	R290 / GWP=3
warranty	24 months
manufacturer	POL-EKO®

all the above technical data refer to standard units (without optional accessories)

* - fluctuation measured in centre of the chamber; in space, variation (K) calculated for chamber as:

$K = \pm (T \text{ average max.} - T \text{ average min.}) / 2$

** - other power supplies on request


1 - depth doesn't include 50 mm of power cable

2 - on uniformly loaded surface

3 - reinforced shelf


4 - reinforced version

OPTIONS AND ACCESSORIES

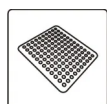

Order number: */PP

Perforated shelf

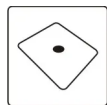

Order number: */PO

Full shelf with hole



Order number: */PW

Reinforced shelf



Order number: KUW GN*/*

Stainless steel cuvettes




Order number: */W

Reinforced version



Order number: OWW/OWW LED

Interior lighting

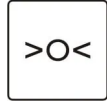

Order number: BRT/*/L or IQ/OQ/PQ

Calibration and IQ, OQ, PQ qualification



Order number: */3.2

Over temperature protection 3.2 class according to DIN 12880



Order number: BPP 12

Battery backup for display



Order number: PORT ALARM


Dry alarm contact



POL-EKO®

(+48) 32 453 91 70
info@pol-eko.com.pl
www.pol-eko.com.pl

POL-EKO A.Polok-Kowalska sp.k.
44-300 Wodzisław Śląski
ul. Kokoszycka 172 C

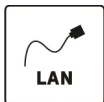

Order number: LANK

LAN cable